


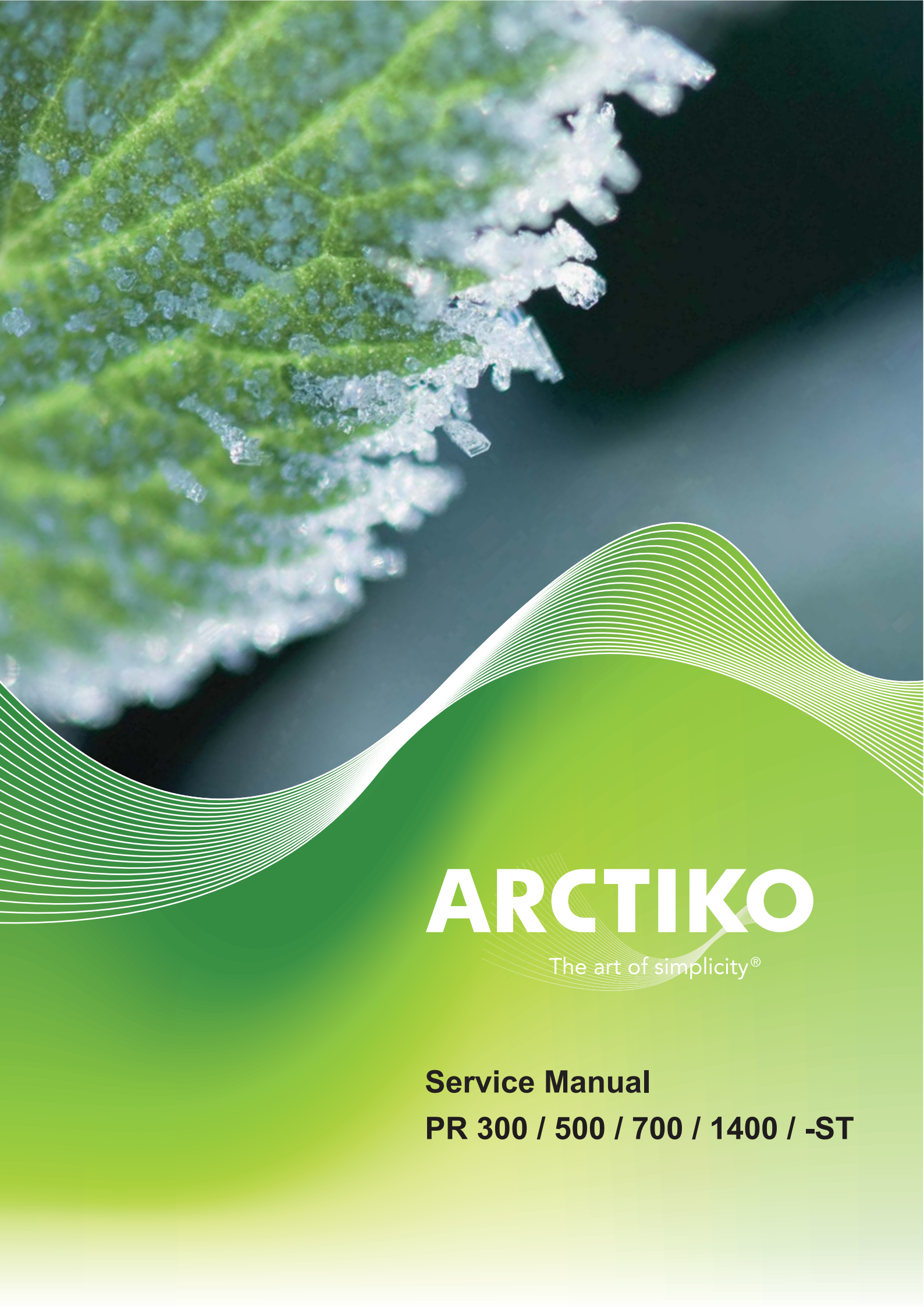

ARCTIKO

The art of simplicity®

Service Manual

PR 300 / 500 / 700 / 1400 / -ST

Index

1. Working safety

- General safety
- Elektrical safety
- Chemical and biological safety

2. Information about the unit

- 2a. General information about the unit
- 2b. Controller G-214
- 2c. Model information

3. Troubleshooting / Repair

- 3a. General operating problems
- 3b. Troubleshooting schedule
- 3c. Opening of the cooling system
- 3d. Evacuation of the refrigerant
- 3e. Changing the dryfilter
- 3f. Changing the compressor
- 3g. Pressure testing and leak detection
- 3h. Vacuum suction / Vacuum check
- 3i. Refrigerant filling with R134a (PR300/500/700/-ST) and R404a (PR1400/-ST)
- 3j. Electrical troubleshooting
- 3k. Electrical fault - compressor
- 3l. HST-start relay, Start capacitor

4. Maintenance

- 4a. Rubber gasket of the door
- 4b. Heated door frame
- 4c. Evaporator, Fan, Heat.element for aut. defrost and heating spiral in the defrost water drain
- 4d. Condenser - fan
- 4e. Door switch
- 4f. Sensor for cold storage temperature
- 4g. Controller - demounting the display
- 4h. Automatical / manual defrost, cleaning and maintenance
- 4i. Noise
- 4j. Battery backup

5. Wiring diagrams and safety data sheet

- 5a. Wiring diagram PR300/500/700/-ST and PR1400/-ST
- 5b. Spare parts list PR300/500/700/1400/-ST
- 5c. Safety data sheet R134a and R404a
- 5d. Controller G-214
- 5e. Default settings

1. Working safety

This page provides general safety information for servicing the unit.

The instruction manual includes additional information about security and cleaning of unit.

1.1 General safety

To prevent personal and equipment damage, the electrical safety / legislation of your country must be observed.

- 1. Never expose the moving parts of the plant to physical overload.
- 2. When repairing plant parts, safety concerns must be respected.
- 3. When using tools, spare parts and equipment, the existing security concerns must be met.
- 4. Maintenance procedures described in this manual must only be performed by trained personnel.

1.2 Electrical safety

Maintenance procedures may imply a voltage hazard. Follow existing safety instructions, including the low voltage and current EU directives and other directives of your country.

Always disconnect the power to the unit before repair or parts replacement is initiated.

Use only electrical equipment which is designed for the unit.

1.3 Chemical and biological safety

- Always observe the requirements for safe handling, as described below.
- Read the instruction manual before using, cleaning and servicing the unit.
- Wear eye protection and gloves when working with refrigerants.
- Skin contact with liquid refrigerants may cause frostbite.
- Good ventilation and air extraction at the work area is required.
- Stay in environments with high concentrations of refrigerant vapors is harmful.
- Low concentrations may cause narcotic effects.
- Symptoms of this may include: headache, shortness of breath, nausea and difficulties in concentration. Possibly: Mobility / consciousness.
- Read the safety data sheet section 5d.

2. Information about the unit

2a. General information about the unit


PR 300 / 500 / 700 / -ST



PR 1400 / -ST



Label plate located on the device

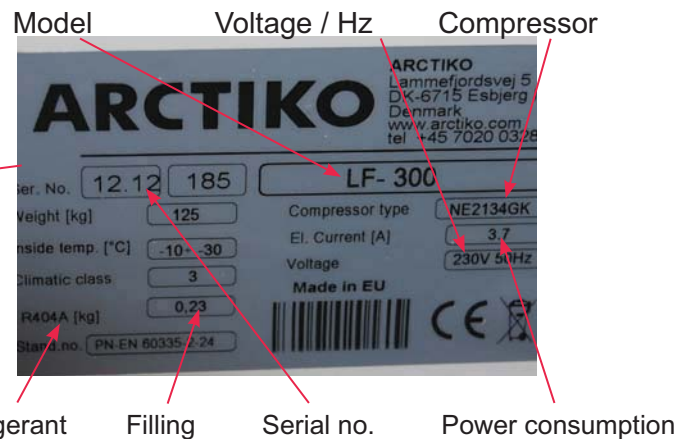

When contacting the supplier it is important to inform the model and serial no.

The label plate is located behind the front panel in the engine compartment.

For technical problems: See troubleshooting schedule (3b).

Furthermore, the device is equipped with heated door frame.

The unit must be defrosted manually / automatically. See section 4h.

2b. Controller G-214

The PCB controller board is located in the engine compartment behind the front panel.

The display is mounted in the front panel.

Data controller:

- 72 hour battery backup for alarms.
- Microprocessor with digital display.
- Temperature display and data acquisition during power failures.
- Visual and audible alarm.
- Adjustable high / low temperature / alarm.
- Integrated data logger.
- Automatic cycle of sensor failure.
- Measurement of the ambient air temperature.
- Direct download / upload from the USB memory stick.
- Display of battery voltage and ambient temperature.
- Integrated memory and data logger.
- RS 485/232 interface.
- Contact for remote alarm.
- Prepared for GSM alarm.

Alarm for:

- High / low temp.
- Voltage fault
- Sensor error
- Instrument error
- Front door

Display G-214



EI-box

PCB Board

Data sheet for controller G-214.
See section 5d.

2c. Model information

Model	PR300 / -ST	PR500 / -ST	PR700 / -ST	PR1400 / -ST
Item code	DAI 0885 / 0886	DAI 0850 / 0851	DAI 0860 / 0861	DAI 0870 / 0871
Temp. range (°C)	+1/+10	+1/+10	+1/+10	+1/+10
Max. ambient temp. (°C)	43	32	32	32
Power supply (V)	230/110*	230/110*	230/110*	230/110*
Frequency (Hz)	50/60*	50/60*	50/60*	50/60*
Fuse (A)	10	10	10	10
AMP (A)	3.29	1.2	1.8	3.2
Power consumption (KWh/24h)	2.4	3.53	5.7	8.13
Compressor	See label plate	See label plate	See label plate	See label plate
No. of compressors	1	1	1	1

Cabinet material on ST-models: Stainless steel

* Special produced units (not standard).

Compressor:

The compressor type can be found on the rating plate in the engine compartment.

The units may be provided with compressors of different brands, but with the same specific properties.

3. Troubleshooting / Repair

3a. General operating problems

Malfunction of the device, which can arise from improper handling of the device:

Moving or long periods of inactivity of the unit

When moving or after long periods of inactivity, you must wait 1-2 hours before the unit is started.

Loading the unit

When loading large quantities of items into the unit, the temperature will rise temporarily.

Repeated door openings

In case of repeated door openings, leaks or faulty gaskets, the humidity will cause rime and possibility of ice formation inside the device and around the door / interior door.

Component failure occurred during operation

Use the troubleshooting schedule (section 3b).

3b. Troubleshooting schedule

Problem	Cause	Recommended equipment	Troubleshooting	Action
The unit does not freeze and the compressor does not start.	The socket or fuse is defect.	Universal Instrument test lamp.	Check the installation.	New fuse / plug.
	Power cable is defect.		Check power cables and wiring in the engine compartment.	Repair errors. Consider using wiring diagram section 5.
	The controller is defect.		Test sensor.	Change controller.
	The sensor is defect.		Test compressor windings.	See section 4l.
	Faulty windings in compressor		Test start relay.	See section 3k.
	Start relay defect.		Use wire-diagram.	See section 3l.
	Fault in el-installation.			See section 5.
The unit does not freeze. The compressor tries to start, but fails.	Low voltage.	Voltmeter	Check voltage.	Find the cause.
	Defect or faulty start relay/capacitor.	See section 3m.	See section 3m.	Change start relay/capacitor.
	Defect compressor windings.	Ohmmeter	Check for electrical failure in the compressor.	Change compr. See section 3f.
	Defect rotor (compressor failure).			Change compr. See section 3f.
	High pressure in the coolin system.	Tapping valve Manometer Filter	Open the system and rinse the system with nitrogen.	Remove 5 cm of the capillary tube. Change filter. See section 3e.
	The unit has been transported horizontally.			The unit must rest a few hours before start.
	Ambient temp. low. Compr.oil too cold.	Eksternal thermometer	Check the temperature.	Find a better placement.
	May be general, if the ventilation is not sufficient.	Eksternal thermometer.	Check, that the compressor not exceeds 70°C.	Find a better placement or ventilation
	Ambient temp. high.	Thermometer	Check ambient temp.	Find a better placement or ventilation

3b. Troubleshooting schedule

Problem	Cause	Recommended equipment	Troubleshooting	Action
The compressor runs constantly. The unit freezes normally or colder.	The sensor is incorrectly installed.			
	Sensor is defect.			
	Controller is damaged.			
	Icing on the frame.		Check the door for leaks.	See section 4a.
	The temperature is set too low.		Check the settings in the controller.	Set new values. See section 5e.
The compressor starts normally, but stops again.	Too high voltage.	Voltmeter	Measure the voltage.	Inform.
	Too high ambient temperature.	Eksternal thermometer	Check the temperature and ventilation.	Better placement or ventilation.
The unit cools too much, too little or not at all. The compressor may run continuously	Leakage in the system.	Elektronic leak-seeking, soapy water, leak seeking spray.	First check the solder points, compr., capacitor etc.	See section 3.
	The fan is not running.	Voltmeter	Check for errors.	Replace the fan. See section 4c.
	The evaporator is blocked by ice.			
The unit cools too much or not at all. Compressor running constant.	Little or no compressor capacity.	Service valve Manometer	Mount the valve and check pressure.	Change compressor. See section 3f.
Insufficient cooling.	The device has recently been filled up with a large amount of heat.			Inform the user.
	In addition, temperature of the air is too high.	Eksternal thermometer	Check temperature.	Find better ventilation or placement.
	In addition, temperature of the air is too low.	Eksternal thermometer	Check temperature.	Find better ventilation or placement.



Examples of troubleshooting with a manometer connected to the process piping of the cooling system.

The system is blocked:

Suction pressure is very low.

Pressure equalization takes place slowly or not at all.

Reason: Icing or dirt in the capillary tube or filter, or blockings in the system due to other reasons.

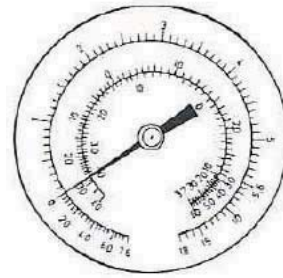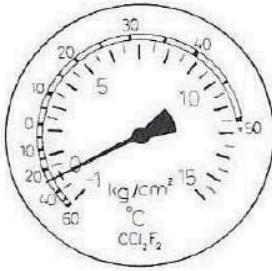


Reason: No compressor capacity.

The suction pressure is too high.

The suction pressure does not change noticeably when the compressor stops.

This can be tested with a volumetric meter.

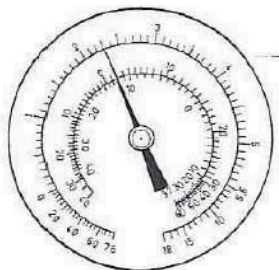

Leakage in the system.


The pressure gauge indicating insufficient suction pressure.

On the pressure side the temperature increases.

Pressure equalization time is short.

countervailing pressure is lower than expected.

Be sure that the pressure is sufficient for a pressure test.

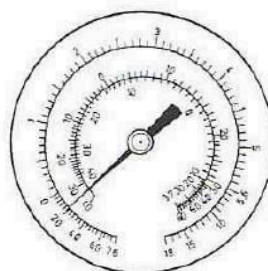

3c. Opening of the cooling system


- To prevent moisture from penetrating under repair, the system must not be open for more than 15 minutes.
- Prepare spare parts components before the system is opened.
- The new connector is soldered to a tube with valve.



- Verify that the new filter packaging is intact and that the filter is not damaged.
- If the packaging has been broken for a longer period, the filter will be moist.

3d. Evacuation of the refrigerant

- It is recommended that the compressor is running at least 25 min. before refrigerant evacuation is undertaken. This ensures that the system is hot, so that the refrigerant easier can be pulled out of the compressor oil.
- A drilling tongs with valve is placed on the process pipe.
- The refrigerant is evacuated.

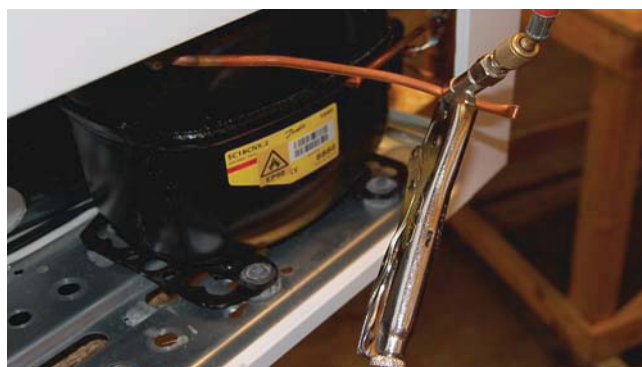

- Install a valve on the pressure tube and process pipe.

Process pipe and pressure tube fitted with a valve.




- The cooling system is evacuated with a two-stage vacuum pump for at least 30 minutes. Both from the pressure piping and process piping.
- Blow the system with nitrogen, both from the suction and the pressure side.

3e. Changing the dry filter

- The dry filter must pick up moisture in the system during operation and function as the capillary tube in front of the capillary tube, in order to avoid clogging. Clogged filter or capillary tube causes a pressure drop, by this the filter becomes colder and ice will block the system. The compressor is overloaded with a possibility of crash.
- The dry filter is placed after the condenser outlet, immediately before the capillary. Any opening of the cooling system, dry filter change to ensure the function and lifetime of the unit.
- It is therefore recommended always to order a new filter drier, along with a new service cylinder with cooling agent..
- To prevent moisture from penetrating during the repair, the drying filter must not be open for more than 10 min.
- After evacuation and blowing with nitrogen, the tubes must be cleaned on both sides of the filter drier, before cutting. By doing this the pipe ends are free from impurities and metal pieces and ready for soldering.
- Cut off the pressure pipe close to the filter with a pipe cutter.




- Cut off the capillary tube with a capillary tube scissors.
- It is recommended to use a special tongs to create a wave on the capillary tube which ensures the correct positioning in the filter.

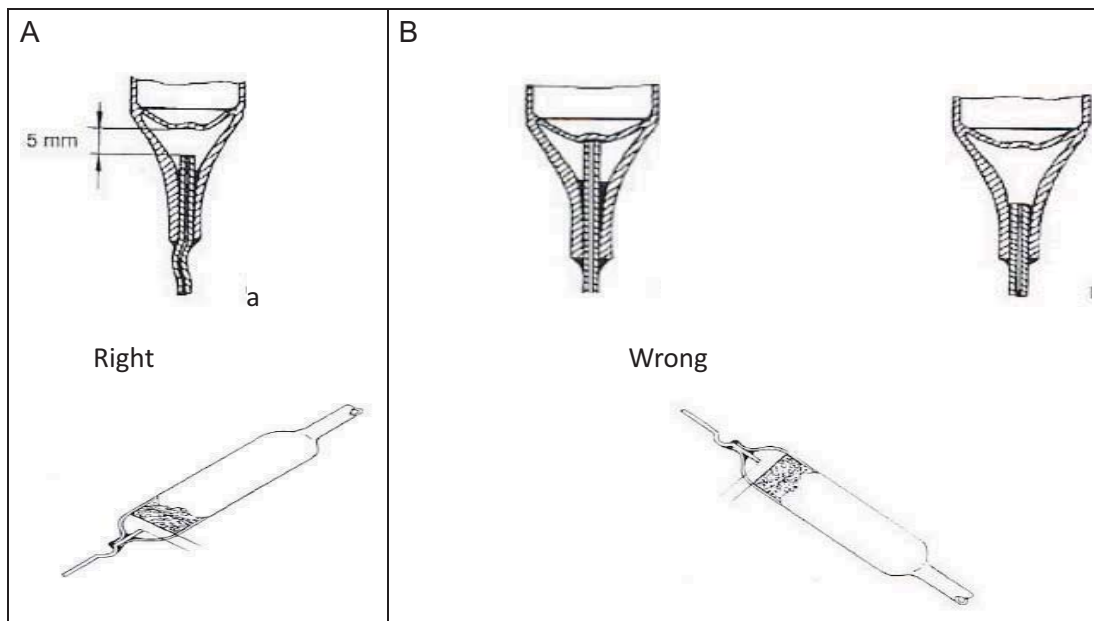

In the figure A the capillary tube is correctly installed in a dry filter. In order to ensure efficient use of the filter, it must be positioned with an inclination of at least 150 and the capillary tube must be underneath.

In the figure B the capillary tube is too close too the web of the filter. It provides great resistance and filter blocking after a short time. Or too far away, so it will be filled with flux or solder.

- Solder the new filter.
- Clean the solder points with a wire brush.

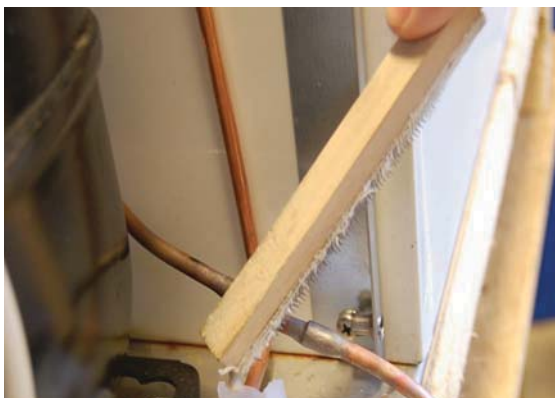

- Flush the system with nitrogen.
- If only the dry filter is to be replaced, the system must be pressure tested. See section 3g.

3f. Changing the compressor

- Evacuate the refrigerant of the system. See section 3d.
- All tubes are cut with a pipe cutter. Do NOT solder of.
- When replacing the compressor the dry filter must also be replaced. See section 3e.
- The suction and pressure pipes of the compressor must be cleaned and cut with a pipe cutter just above the soldering pipe.
- Demount the compressor.



- Install the new compressor, which is prepared as described before. (3c)
- In order to avoid decomposition of refrigerants which may remain in the system during the soldering operation, the system is blown with dry nitrogen, respectively, from the suction side respectively pressure side.
- Solder the tubes on to the compressor.
- Pressure test the system.

3g. Pressure testing and leak detection

- Pressure test the system with nitrogen. The pressure must be 10 bar during a period of 2 hours.
- Check the solder points with soapy water or leak detection spray.
- In case of any leaks, empty the system for nitrogen.
- Repair and pressure test again.
- During operation leak detection must be performed on the pressure side with the compressor running. The suction side must be tested while the compressor is stopped, and the pressures are equalized.
- The nitrogen must slowly be emptied from the system.

3h. Vacuum suction / Vacuum check

A vacuum is made on two locations on the system, one on the process tube and one on the pipe branch used for testing located on the pressure tube with an explosion-proof 2-stage vacuum pump.



- The pressure should be 2 mBar or lower for at least 12 hours.
- The system is now ready for the filling of the refrigerant.

3i. Refrigerant filling with R134a / R404a

Before filling the refrigerant, a vacuum must be sucked. See section 3h.

Note that R134a / R404a is to be supplied from a bottle containing a specific amount of refrigerant.

Common procedures for other refrigerants can not be used.

Further it is important that the entire contents of the bottle is introduced into system to ensure optimal performance.

The contents of the bottle must be mixed before use by shaking the bottle for at least 1 min. before the filling.

In order to perform the filling operation properly, it is necessary to use a capillary tube with the valve before the filling site, instead of a conventional feeding tube. The reason is that the capillary tube ensures that the entire quantity of refrigerant is introduced to the system, If an ordinary filler is used, it must be as short as possible.



The refrigerant cylinder is weighed before use.

Cylinder with manometer



Vacuum suction of the filling hose



Loading procedure

Make vacuum. See section 3h.

The cylinder with manometer is connected to the coupler of the process piping.

Remember first to make vacuum in the filling hose to empty it.

Fill the refrigerant into the system until the pressure gauge on the pressure tube shows approx. 5 bar. Start the compressor and let it run for 2-3 min. Open to the cylinder again to get out the remaining refrigerant.

When the pressure in the bottle shows 1 bar, the bottle is empty.

It is important that the valve is closed immediately, in order to avoid that the refrigerant runs back into the cylinder.

Now the cooling system is filled with refrigerant and only needs to be closed.

Disconnect the filling hose with pressure gauge from the process piping.

Push the process piping and pressure pipes together with a special tongs.

Cut off valve with a pipe cutter.



Closing the process tube on systems with flammable refrigerants must be performed by professionals with a LokRing end sleeve. It is recommended NOT to solder system.

Clean the pipe with steel wool or emery cloth. Use rotational movement in order to avoid scratching the length of the tube.

Apply LokPreb to the end of the tube.

Turn the end sleeve 360°, so that the LokPreb is placed around the tube.

Squeeze the end sleeve together with the special tool as shown in the illustration..

Start the unit and control the function of the device.


3j. Electrical troubleshooting

Before starting systematic troubleshooting, check the following:


- Proper voltage is present.
- Electrical equipment used, is suitable for the compressor.
- The wires are properly installed according to wiring diagram.
- There is no transition between the live parts and ground.

Control of the transition to ground.

In order to check whether there is a transition between the conducting components and the chassis, the insulation resistance is tested by using an Ohmmeter or the second high-voltage equipment. An Ohmmeter can be used for a rough check. The Ohmmeters one clamp is placed on the compressors ground terminal. The second clamp on the compressor shared leg (see sketch). Is the Ohmmeter showing a transition, the compressor is replaced.

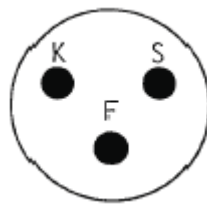

Plug on the compressor

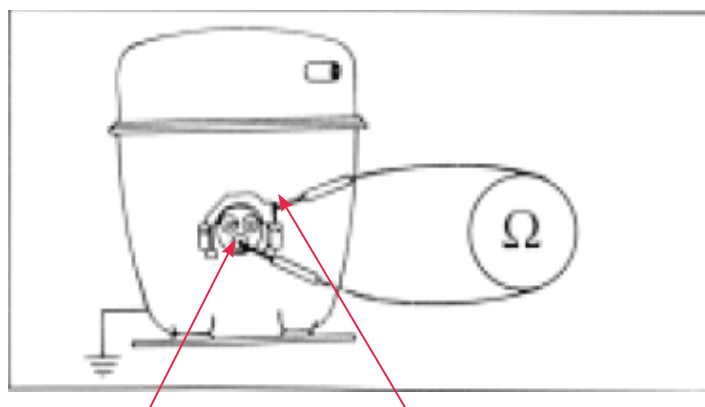


Shared leg

Ground terminal

Uses high voltage test systems with flammable refrigerants, ensure that the system is completely drained of coolant. as this can cause sparks during the test.


3k. Electrical fault - compressor

Remove the start relay from the compressor and use an Ohmmeter to test the compressor head and start recovery. Connect the Ohmmeter between the common connector and start winding, then the main winding to debug. See the picture below.

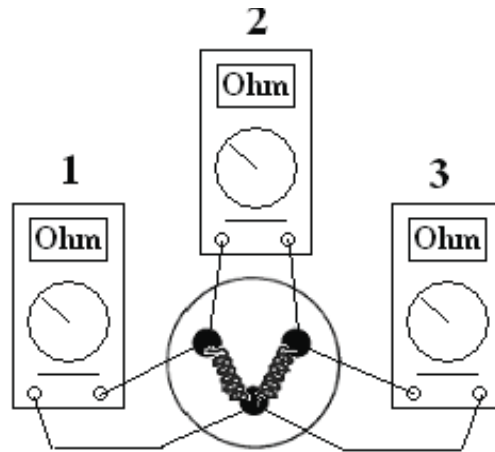

Dimensions of the compressor connections to determine whether driving and start winding are intact.

This is done by measuring the resistance of these. Carry out three resistance measurements on compressor penetrations. The sum of the two smallest measurements, the intact windings will be equal to the highest. Find if the current values of the compressor data sheet.

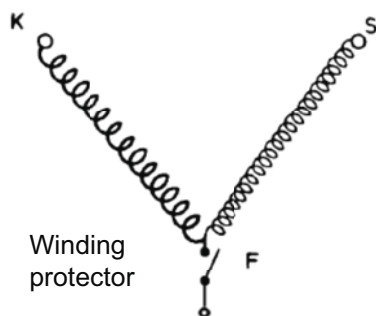

Ohm-measurement

By measuring the compressors with built winding protection it is extremely important to measure between K and S, as this measure will tell if this is burnt windings or an open winding protector.


Since it is relatively small value to be measured, it would be a great advantage to use a digital Ohm-meter.


If the measurement indicates error, the compressor is replaced.

3l. HST-start relay, start capacitor


Connect a voltmeter between terminals 10 and 13 at the start of the relay. If no voltage is present, there is either a faulty cable or also the start relay is defect.

Connect a voltmeter between terminals 10 and 14. If no voltage is present, the temperature sensor or its cable is defect.

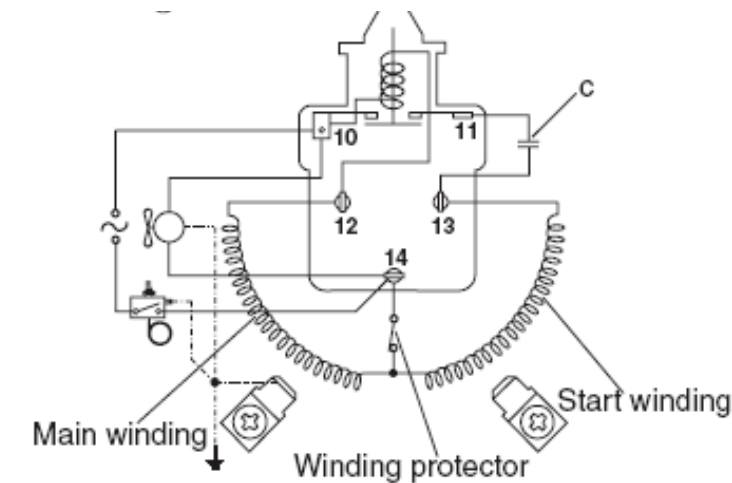

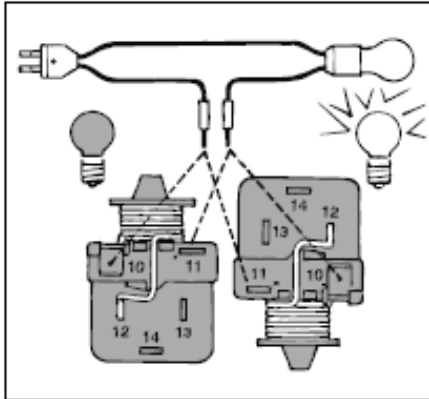

A relay can be checked with an indicator light as shown in the sketch maintain end. The relay is usually ok, when in standing position is not lit, and when it is lit when the relay is upside down.

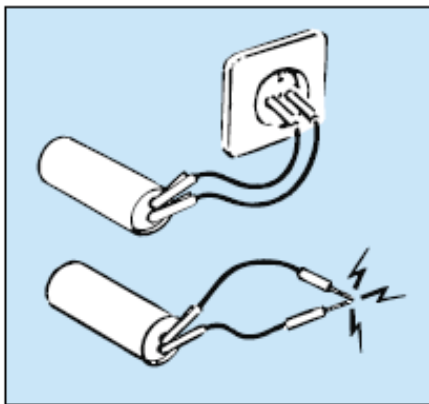


The start capacitor can be checked by injecting the normal voltage for a few seconds and than shorting it. When sparks are generated, the capacitor is in order.

When replacing the compressor the start relay and the start capacitor should also be replaced.

The start relay is placed and clamped on the three terminal pins on the compressor.

The voltage is connected to terminals 10 and 14.

The start capacitor is connected and clamped to terminals 11 and 13.

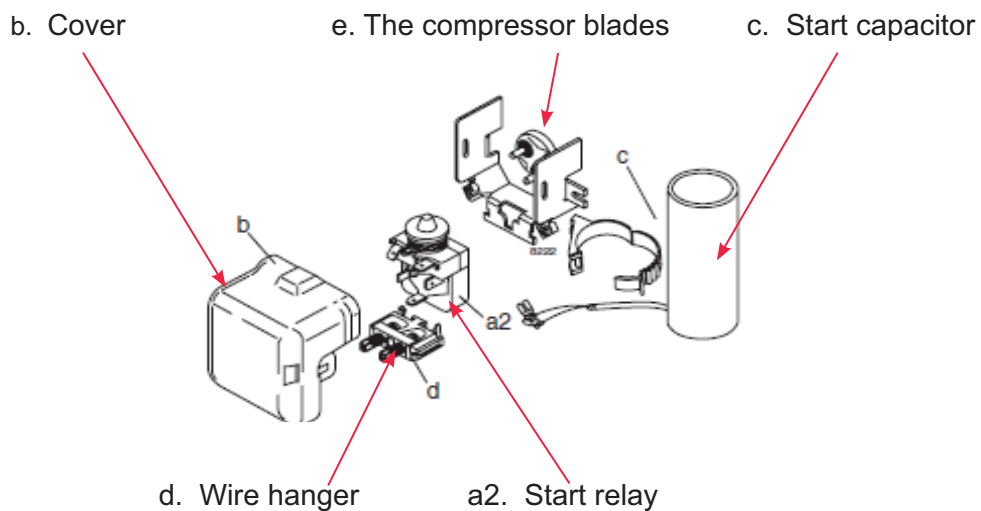

4. Maintenance

4a. Rubber gasket of the door

Rub the gasket with a suitable product, so that it stays flexible and the optimal function and tightness is secured.

Make sure the door closes tightly so that no ice forms around the door frame.

- Open the door and place a piece of paper. Close the door.
- Pull the paper.
- If the paper can not be pushed back, the adjustment is OK.
- The test is made at both the left and the right side.
- In case that the gasket do not close properly, see section 4b.



The door can be adjusted in the slots in the top hinge of the door in the motor compartment.



Adjusting

4b. Heated door frame

The heaters are installed in the door frame .

Electrically, the cable (220V) is supplied by the PCB board and assembled in the controller box.
See the wiring diagram section 5a.

The cable is introduced to the door frame on the front of the engine compartment behind the front panel.



4c. Evaporator, Fan, Heat.element for aut. defrost and heating spiral in the defrost water drain

The evaporator is placed at the top of the cold storage.

The fan runs continuously when the unit is connected and ensure good air distribution in the cold storage.

The fan stops when the door is opened and during the automatic defrost.


To avoid ice formation in the defrost water tray, a heater is mounted in the drain pipe to the drain tray.



Fan LR 500 - 700 - 1400 (2 pcs.)



Fan LR 300

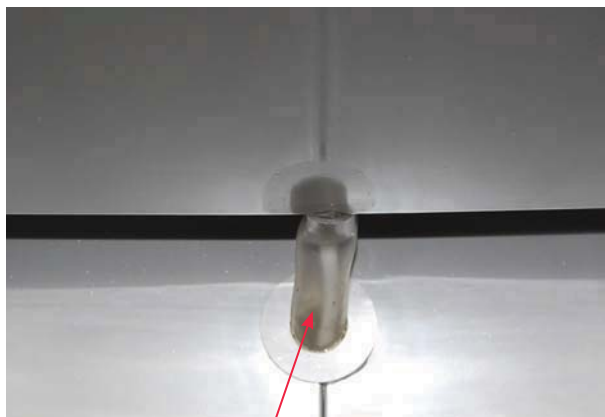


Heater in the drain pipe


To get into the evaporator, the sensor for evaporator temperature and the heating element the drain is removed. Then remove the drain tray (4 screws).



Remove the drain pipe

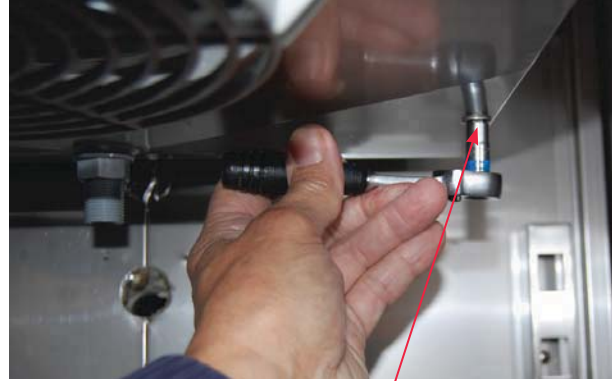


Remove the drain tray

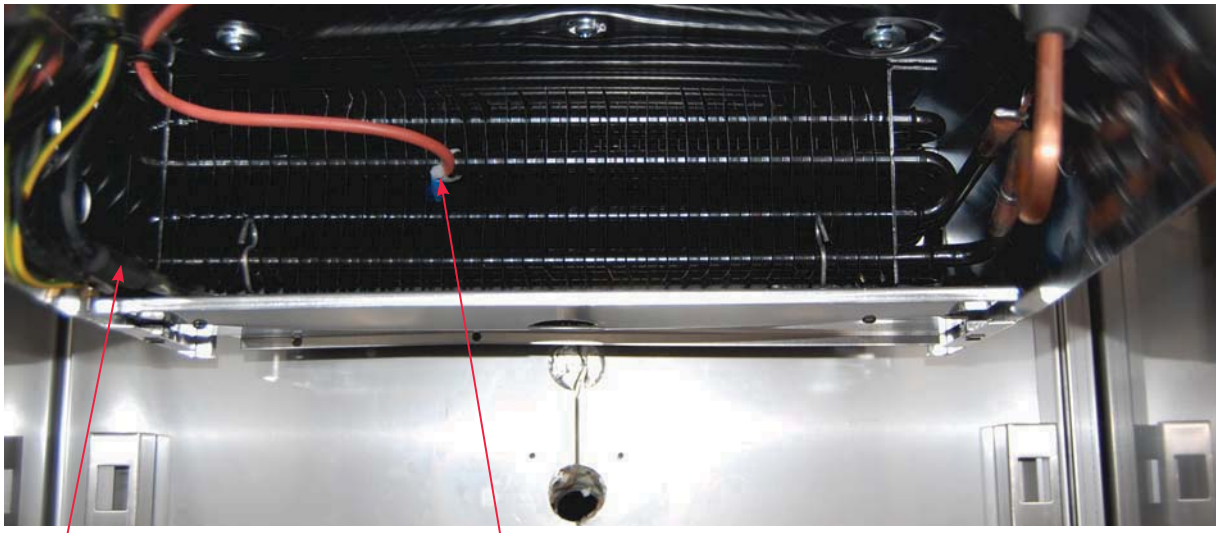

Heating element

The sensor for evaporator

Sensor lead, wire for ventilation and heater as seen from the top of the motor compartment.

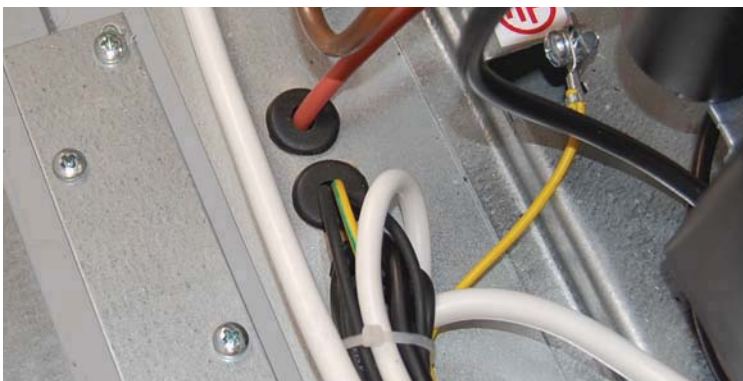

4d. Condenser - Fan

The condenser is placed in the engine department of the uit behind the front panel.

Check that the bolts are well tightened and that the rubber gaskets at the fan are not hard or defective. It is important that the cooling ribs are kept clean and the ambient temperature does not exceed 25°C.

Physically, the fan is placed on the frame of the capacitor.


Electrically, the fan is connected to terminal 10 and 14 on the start relay. See wiring diagram section 5a.

If the ventilator stops, it can cause insufficient cooling of condenser and compressors.

Higher condensing temperature means that the compressor must work more with less efficiency as a result. At worst, the compressors stops due to high temperature.




Fan

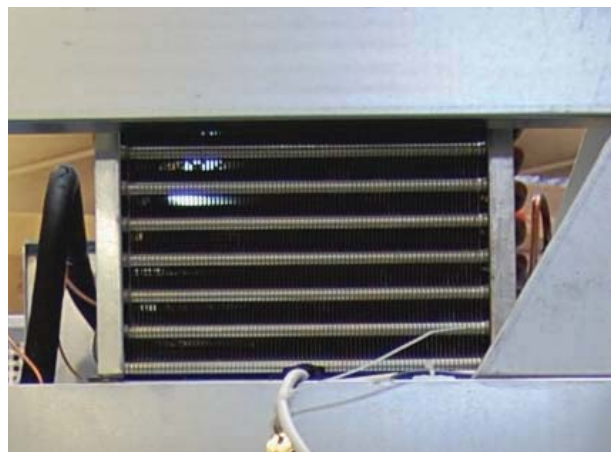


Cooling ribs


4e. Door switch

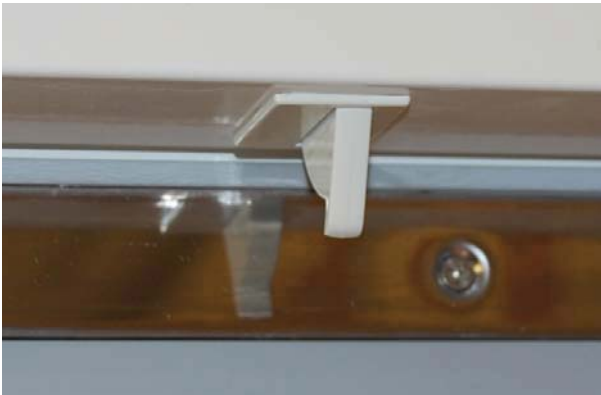

Door switch under the top panel

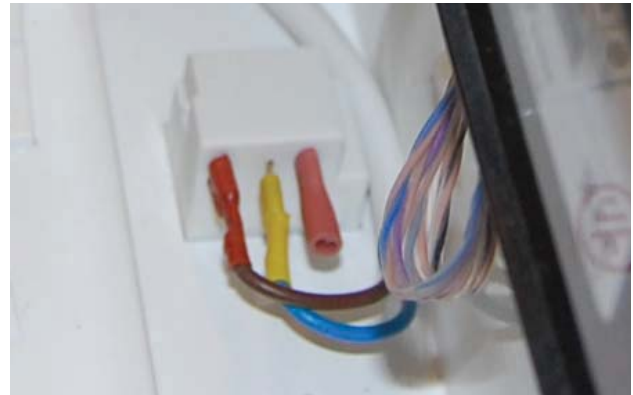

Door switch in the motor compartment

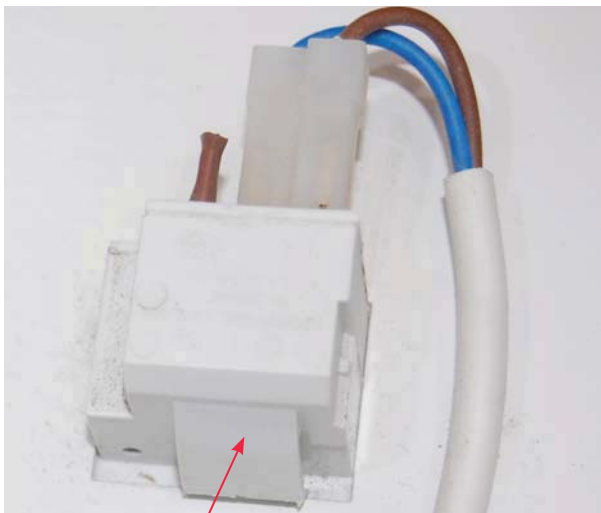

Push here when remounting

The door switch function controls:

Start/stop of the evaporator fan.

On/off light.

Alarm for open door.

The door switch is clipped under the door panel.

When replacing remove wire.

The flaps are pushed together and the switch is pressed out.

A new switch is mounted in reverse order.

4f. Sensor for cold storage temperature

Sensor for evaporator tube temperature:

The sensor is the type of PT 1000.


The sensor must be introduced from the top of the motor compartment as the sensor for the room temperature.

The sensor is placed at the top of the evaporator in the cold storage covered by the drain cover.

To remove the sensor, remove the cover. (See section 4c)

The sensor is fastened with strips.

The sensor is connected to the PCB board input. (See wire-diagram section 5a)

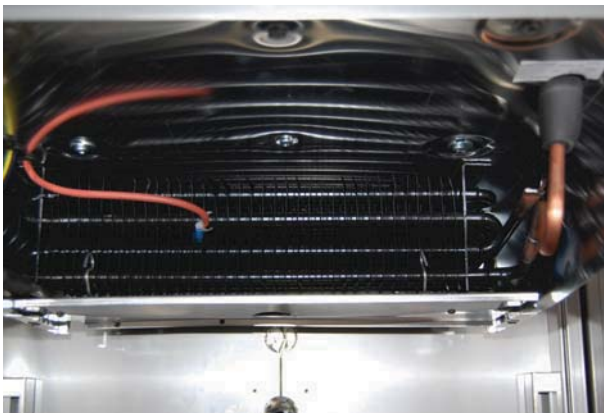

Evaporator



Sensor wire

Sensor for cold storage temperature:

The PT 1000 sensor for cold storage temperature is placed on the right side fastened with strips.

The sensor must be introduced from the top of the motor compartment.

The holes are sealed with putty.

Electrically, the sensor is connected to the PCB board input (Temp. sensors - Chamber).


See wire-diagram section 5a.

Change of probe:

- Disconnect the power supply.
- Disconnect the electrical connection of the sensor in the CBS modul.
- Remove the strips.
- Remove the sensor track.
- Pull the sensor into the motor compartment.
- Mount the new sensor in reverse order.

Sensor track

PT 1000 sensor

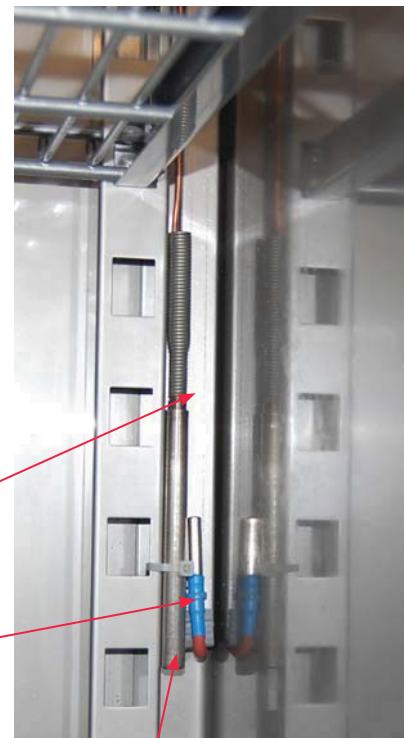

Sensor for possible Chart-recorder

4g. Controller - demounting the display

Demount of display

- Disconnect the wire connector on the bottom of the display.
- Squeeze the mounting devices at the bottom and top together and push out the display.



Mounting device

Wire connectors

Mounting device

4h. Automatical / manual defrost, cleaning and maintenance

There will be some ice formation on the sub lids of the unit.

The ice must be removed with a scraper on a regular basis. (Remember safety gloves)

Ice on the inside of the cabinet must be removed as required.

This is done by manual defrost:

- Remove the contents of the unit.
- Unplug the unit.
- Leave the door open until the ice can be removed.
- Wipe the interior cabinet and shelves with a cloth until it is dry.

Ice formed on the inside of the door, must be removed with a wooden / plastic scraper. This ice is formed due to the low temperatures inside the cabinet, and if the ice is not removed, it will cause that the door does not close tightly.

Automatic defrost cycles in G 214 controller.

As default the controller will run automatic defrost every 6 hours.

1. The automatic defrost cycles start.
2. The compressor stop.
3. The evaporator-fan stop.
4. The evaporator-heat element starts heating up the evaporator.
5. The evaporator-sensor is measuring the temperature in the evaporator.
6. When the temperature in the evaporator reach -5°C., the heat element stop. parameter d9
7. The ice which was build up in the evaporator, will get into liquid form, and start floating out of the chamber, into the water tray on the back of the unit.
8. The parameter c3 is the time for the water can drip of the evaporator when the heat element has switch off.
9. After the C3 dripping time, the compressor start running again.
10. The parameter C13 will determine when the evaporator fan is starting, if the temperature d8 has been reached.
11. When the temperature in the evaporator reach -5°C., the evaporator-fan start running again.
12. The automatic defrost cycles end.

C7 is the paramet which block the display to measure the temperature in the chamber. Examble: If the display show -28° before the defrost starts, the display will show -28°C during the defrost cycles and 5 minutes (c7) after the defrost is finished. After the 5 minutes the display will show the actually temperature again.

The parameter c14 is start blocking time defrost cycles.

The parameter c15 is the end time of blocking defrost cycles.

Examble : If the customer do not want the unit to be defrosted in normally working hour, 8 o'clock until 16 o'clock.

c14 : 8 and c15 :16, now the unit will not do any defrost cycles in between these hours.

During the automatic defrost cycles, the displayed temperature will show the temperature in the air.

This is not the actually temperature which are in the samples of medicine or similar product.

The temperature will be dropping approx. 2°C in the samples, during the automatic defrost.


4i. Noise

The following must be checked in case of noise:

- The compressor is not tightened well.
- The rubber feet are defect or too hard.
- The compressor process tube is too long.
- The fan blades are damaged or dirty.
- The fan is not tightened well.
- The fan bushings are worn.
- The freezer is not positioned on a flat surface.
- Long pipes can cause resonance vibrations.



Rubber feet

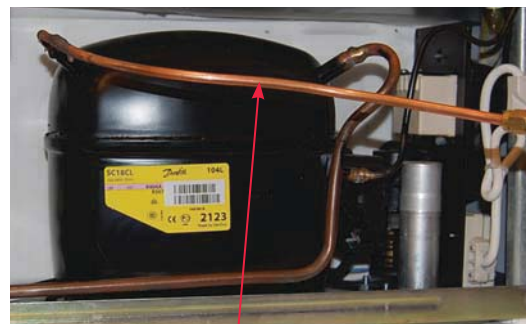

Process tube (too far)

4j. Battery backup

The battery (12V 7AH) is located in the engine compartment behind the front panel.

It has a lifespan of 3 years, but should be checked every six months.


Check the voltage with a voltmeter.

The controller G-214 has built-in measurer for the battery voltage.




Backup-battery


5a. Wiring-diagram PR 300 / 500 / 700 / -ST

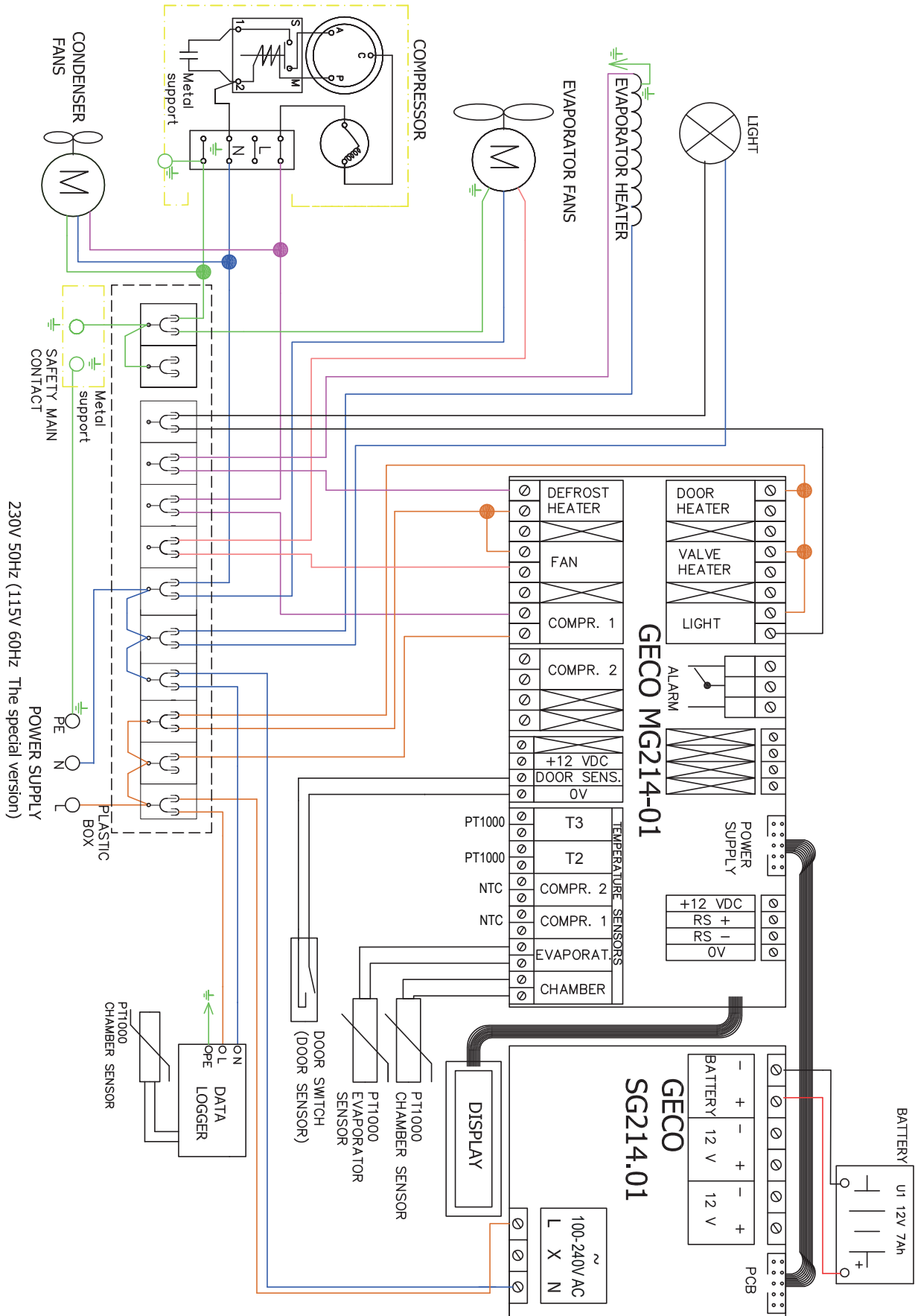

5a. Wiring-diagram PR 1400 / -ST

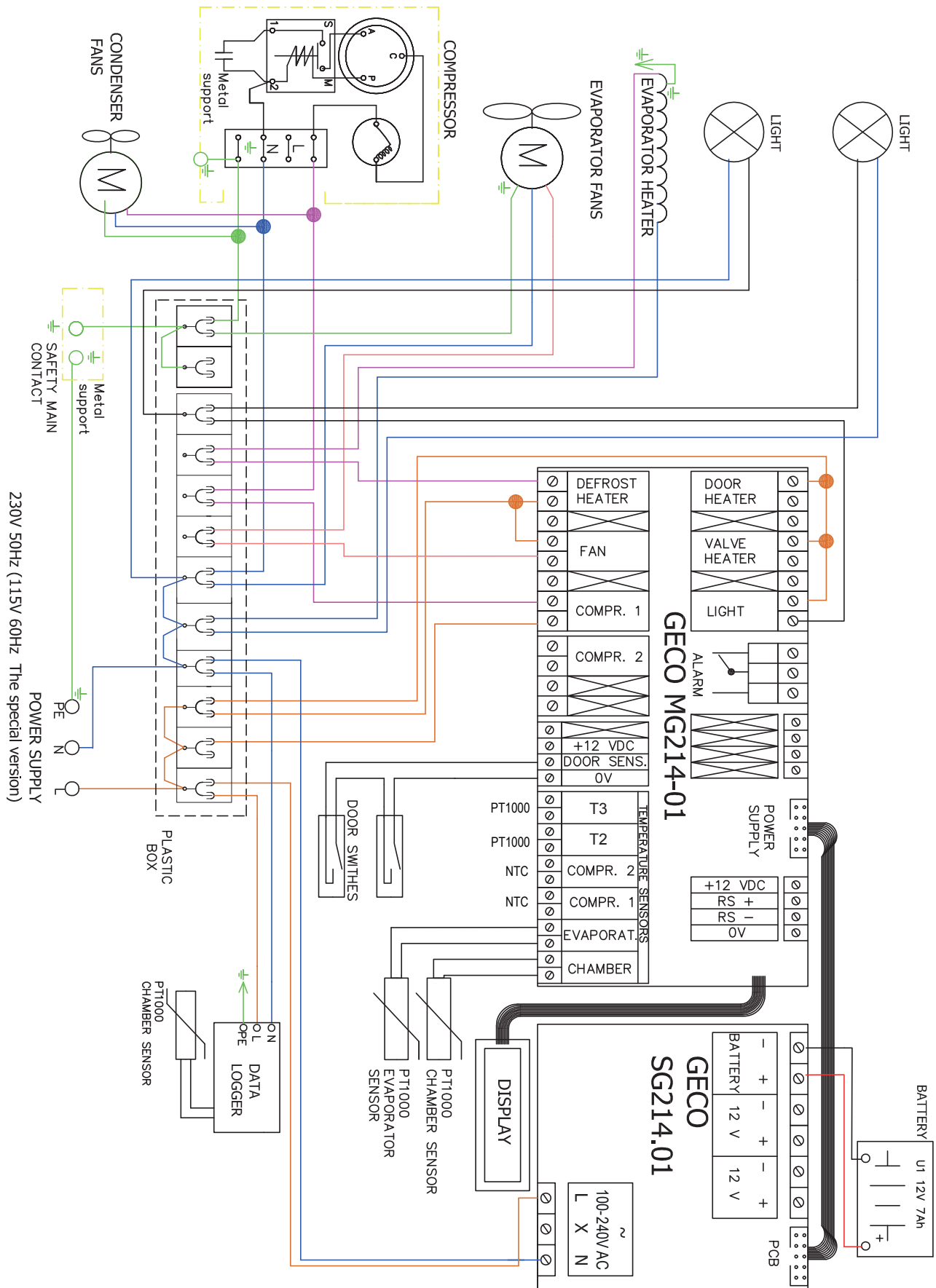


5b. Spare parts list PR 300/500/700/-ST

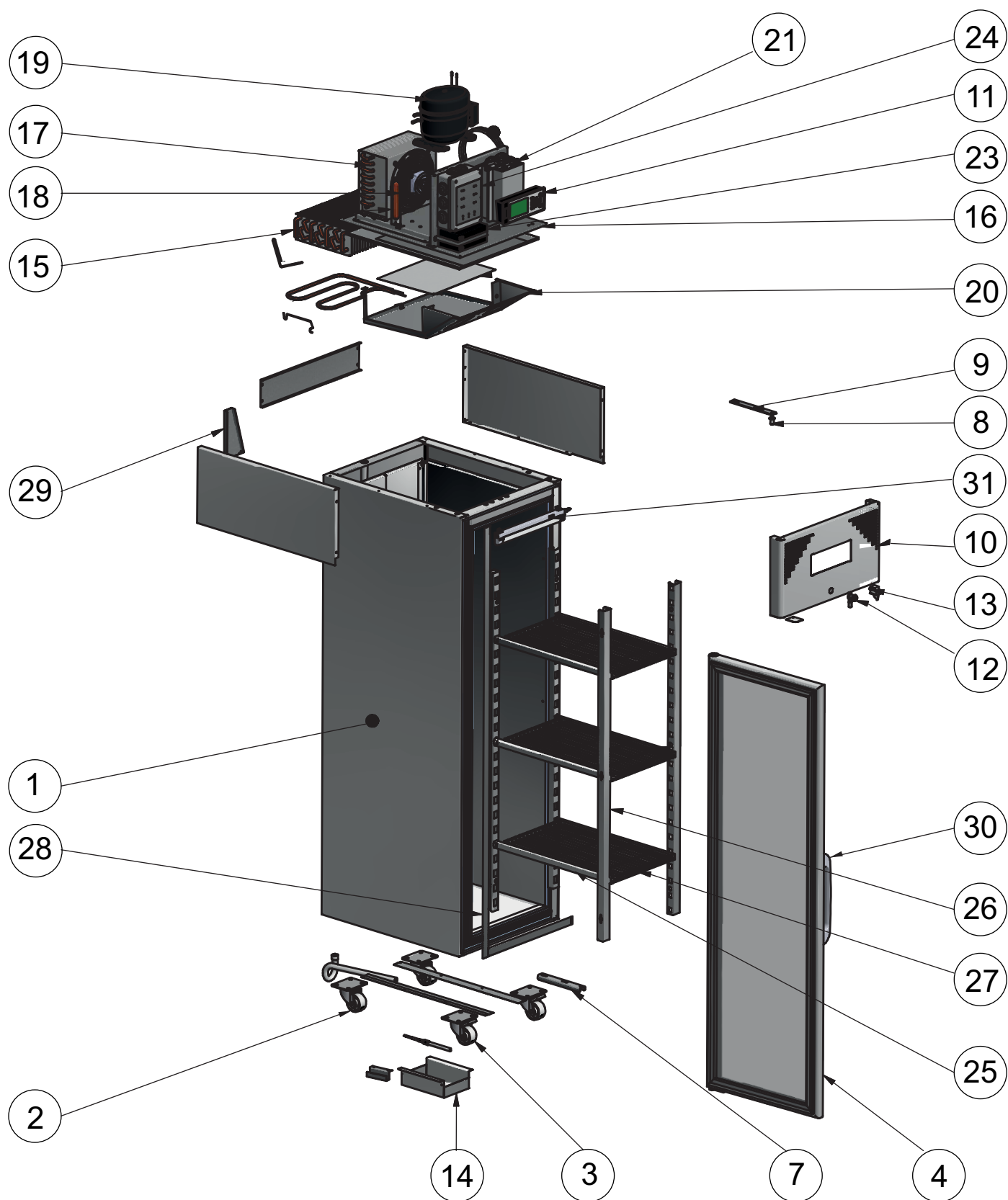

5b. Spare parts list PR 300/500/700/-ST

Pos. no.	Code no.	Part	Recommended
1		Cabinet stainless steel	
1		Cabinet white	
2	5053006-01	Back castor swivel without brake	
3	5053007-01	Front castor swivel with brake	
4		Door	
7	5031034-01	Door lower support	
8	5031012-01	Spring hinge	
9	5031041-01	Door upper support	
10		Front panel Stainless steel	
10		Front panel White	
11	5043007-01	Display G-214	1
11a	5031028-01	Installation plate for display	
12	5051015-01	Lock with key	
13	5019022-01	Fan switch	
14	5031036-01	Defrost water tank	
15		Evaporator	
16	5031138-01	Condenser and compressor support	
17	5071041-01	Condenser	
18	5010007-01	Condenser fan incl. fan blade	
19		Compressor	
20		Evaporator body	
21	5049013-01	Electrical box with support	
23	5044001-01	Battery BP 7-12 (12V, 7AH)	1
24	5041006-01	PCB Board G-214 Controller	1
25	5062009-01	Shelf guide	
26	5062002-01	Shelf support front	
27		Shelf	
28	5062002-01	Shelf support back	
29	5030081-01	Drain	
30	5051005-02	Handle	
31	5019117-01	LED-Light	
-	5041007-01	Power supply G-214	1
-	5019139-01	Power supply support	
-	5042016-01	Probe NTC for compressor	1
-	5042003-01	Probe for chamber (PT1000)	1
-	Finans no.	Refrigerant	
-	5072000-01	Cylinder deposit for refrigerant	
-	5079000-01	Dry filter	
-	5049007-01	Box for PCB board connection.	

5b. Spare parts list PR 1400 /-ST

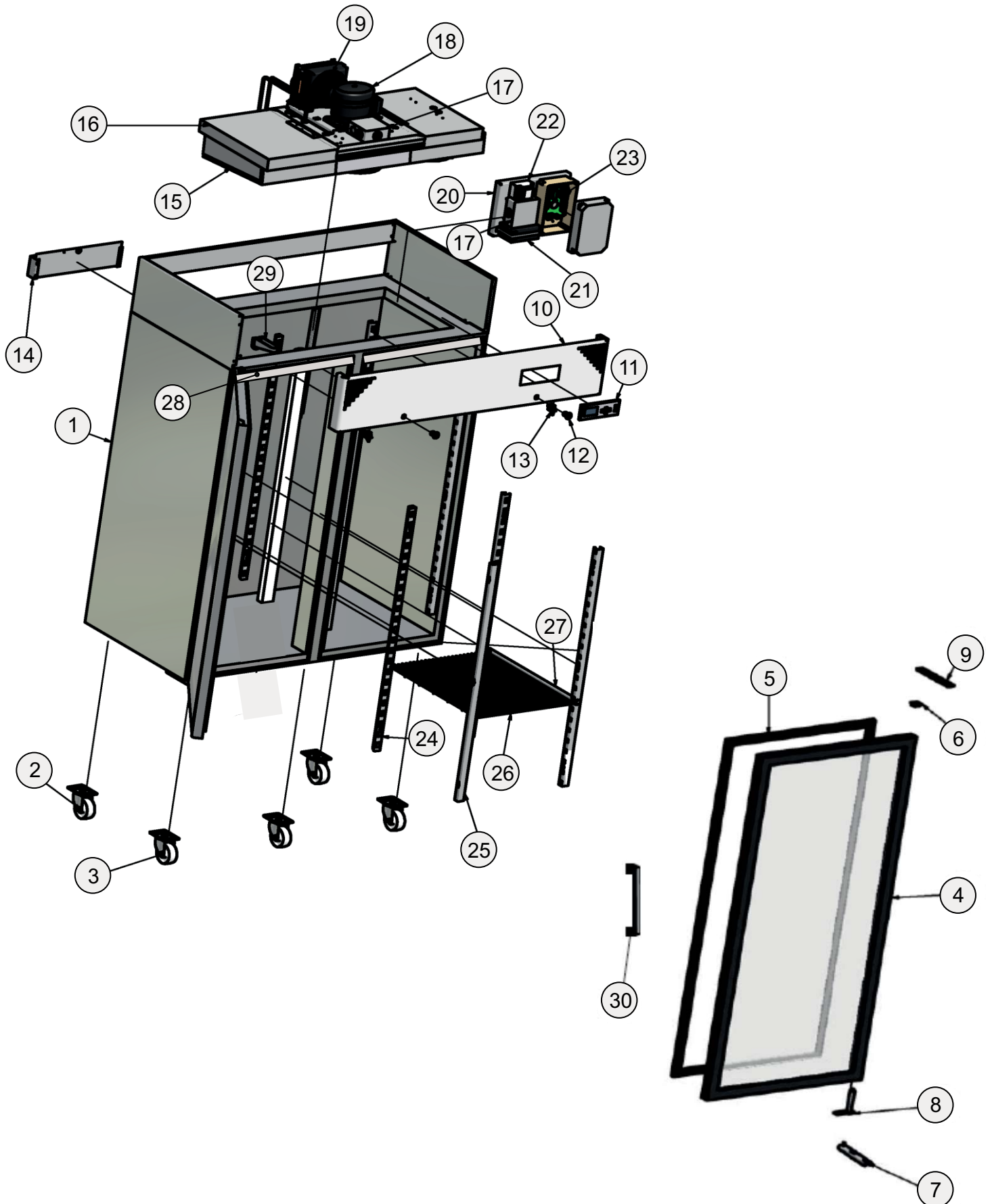

5b. Spare parts list PR 1400 /-ST

Pos. no.	Code no.	Part	Recommended
1		Cabinet stainless steel	
1		Cabinet white	
2	5053006-01	Back castor swivel without brake	
3	5053007-01	Front castor swivel with brake	
4	5052031-01	Door	
5	5052013-01	Gasket	
6	5031033-01	Door sleeve	
7	5031034-01	Door lower support	
8	5031012-01	Spring hinge	
9	5031041-01	Door upper support	
10	5050027-02	Front panel Stainless steel	
10	5050082-01	Front panel White	
11	5043007-01	Display G-214	1
11a	5031028-01	Installation plate for display	
12	5051015-01	Lock with key	
13	5019022-01	Fan switch	
14	5031036-01	Defrost water tank	
15	5073018-01	Evaporator	
16	5031138-01	Condenser and compressor support	
17	5049006-01	Box for 230 V connections	
18		Compressor	
19	5071035-01	Condenser	
19a	5010007-01	Condenser fan incl. fan blade	1
21	5044001-01	Battery BP 7-12 (12V, 7AH)	1
22	5041007-01	Power supply G-214	1
23	5041006-01	PCB Board G-214 Controller	1
24	5062001-01	Sheld support back	
25	5062001-01	Shelf support front	
26	5060002-01	Shelf	
27	5062010-01	Shelf guide	
28	5019119-01	LED-Light	
29	5030081-01	Drain	
30	5051005-02	Handle	
-		Evaporator body	
-	5049013-01	Electrical box with support	
-	5019139-01	Power supply support	
-	5042016-01	Probe NTC for compressor	1
-	5042003-01	Probe for chamber (PT1000)	1
-	Finans no.	Refrigerant	
-	5072000-01	Cylinder deposit for refrigerant	
-	5079000-01	Dry filter	
-	5049006-01	Box for PCB bord connection.	

5c. Safety data sheet for R-134A

Material Safety Data Sheet

R-134A

1. CHEMICAL PRODUCT AND COMPANY IDENTIFICATION

PRODUCT NAME: R-134A
DISTRIBUTOR: National Refrigerants, Inc.
661 Kenyon Avenue
Bridgeton, New Jersey 08302

FOR MORE INFORMATION CALL:
(Monday-Friday, 8:00am-5:00pm)
1-800-262-0012

IN CASE OF EMERGENCY CALL:
CHEMTREC: 1-800-424-9300

2. COMPOSITION / INFORMATION ON INGREDIENTS

<u>INGREDIENT NAME</u>	<u>CAS NUMBER</u>	<u>WEIGHT %</u>
1,1,1,2-Tetrafluoroethane	811-97-2	100

Trace impurities and additional material names not listed above may also appear in Section 15 toward the end of the MSDS. These materials may be listed for local "Right-To-Know" compliance and for other reasons.

3. HAZARDS IDENTIFICATION

EMERGENCY OVERVIEW: Colorless, volatile liquid with ethereal and faint sweetish odor. Non-flammable material. Overexposure may cause dizziness and loss of concentration. At higher levels, CNS depression and cardiac arrhythmia may result from exposure. Vapors displace air and can cause asphyxiation in confined spaces. At higher temperatures, (>250°C), decomposition products may include Hydrofluoric Acid (HF) and carbonyl halides.

POTENTIAL HEALTH HAZARDS

SKIN: Irritation would result from a defatting action on tissue. Liquid contact could cause frostbite.

EYES: Liquid contact can cause severe irritation and frostbite. Mist may irritate.

INHALATION: R-134A is low in acute toxicity in animals. When oxygen levels in air are reduced to 12-14% by displacement, symptoms of asphyxiation, loss of coordination, increased pulse rate and deeper respiration will occur. At high levels, cardiac arrhythmia may occur.

INGESTION: Ingestion is unlikely because of the low boiling point of the material. Should it occur, discomfort in the gastrointestinal tract from rapid evaporation of the material and consequent evolution of gas would result. Some effects of inhalation and skin exposure would be expected.

DELAYED EFFECTS: None Known

Ingredients found on one of the OSHA designated carcinogen lists are listed below.

<u>INGREDIENT NAME</u>	<u>NTP STATUS</u>	<u>IARC STATUS</u>	<u>OSHA LIST</u>
No ingredients listed in this section			

4. FIRST AID MEASURES

SKIN: Promptly flush skin with water until all chemical is removed. If there is evidence of frostbite, bathe (do not rub) with lukewarm (not hot) water. If water is not available, cover with a clean, soft cloth or similar covering. Get medical attention if symptoms persist.

EYES: Immediately flush eyes with large amounts of water for at least 15 minutes (in case of frostbite, water should be lukewarm, not hot) lifting eyelids occasionally to facilitate irrigation. Get medical attention if symptoms persist.

INHALATION: Immediately remove to fresh air. If breathing has stopped, give artificial respiration. Use oxygen as required, provided a qualified operator is available. Get medical attention immediately. DO NOT give epinephrine (adrenaline).

INGESTION: Ingestion is unlikely because of the physical properties and is not expected to be hazardous. DO NOT induce vomiting unless instructed to do so by a physician.

ADVICE TO PHYSICIAN: Because of the possible disturbances of cardiac rhythm, catecholamine drugs, such as epinephrine, should be used with special caution and only in situations of emergency life support. Treatment of overexposure should be directed at the control of symptoms and the clinical conditions.

5. FIRE FIGHTING MEASURES

FLAMMABLE PROPERTIES

FLASH POINT:	Gas, not applicable per DOT regulations
FLASH POINT METHOD:	Not applicable
AUTOIGNITION TEMPERATURE:	>750°C
UPPER FLAME LIMIT (volume % in air):	None*
LOWER FLAME LIMIT (volume % in air):	None*
	*Based on ASHRAE Standard 34 with match ignition
FLAME PROPAGATION RATE (solids):	Not applicable
OSHA FLAMMABILITY CLASS:	Not applicable

EXTINGUISHING MEDIA:

Use any standard agent – choose the one most appropriate for type of surrounding fire (material itself is not flammable)

UNUSUAL FIRE AND EXPLOSION HAZARDS:

R-134A is not flammable at ambient temperatures and atmospheric pressure. However, this material will become combustible when mixed with air under pressure and exposed to strong ignition sources.

Contact with certain reactive metals may result in formation of explosive or exothermic reactions under specific conditions (e.g. very high temperatures and/or appropriate pressures).

SPECIAL FIRE FIGHTING PRECAUTIONS/INSTRUCTIONS:

Firefighters should wear self-contained, NIOSH-approved breathing apparatus for protection against possible toxic decomposition products. Proper eye and skin protection should be provided. Use water spray to keep fire-exposed containers cool.

6. ACCIDENTAL RELEASE MEASURES

IN CASE OF SPILL OR OTHER RELEASE:

(Always wear recommended personal protective equipment.)

Evacuate unprotected personnel. Protected personnel should remove ignition sources and shut off leak, if without risk, and provide ventilation. Unprotected personnel should not return until air has been tested and determined safe, including low-lying areas.

Spills and releases may have to be reported to Federal and/or local authorities. See Section 15 regarding reporting requirements.

7. HANDLING AND STORAGE

NORMAL HANDLING:

(Always wear recommended personal protective equipment.)

Avoid breathing vapors and liquid contact with eyes, skin or clothing. Do not puncture or drop cylinders, expose them to open flame or excessive heat. Use authorized cylinders only. Follow standard safety precautions for handling and use of compressed gas cylinders.

R-134A should not be mixed with air above atmospheric pressure for leak testing or any other purpose.

STORAGE RECOMMENDATIONS:

Store in a cool, well-ventilated area of low fire risk and out of direct sunlight. Protect cylinder and its fittings from physical damage. Storage in subsurface locations should be avoided. Close valve tightly after use and when empty.

8. EXPOSURE CONTROLS / PERSONAL PROTECTION

ENGINEERING CONTROLS:

Provide local ventilation at filling zones and areas where leakage is probable. Mechanical (general) ventilation may be adequate for other operating and storage areas.

PERSONAL PROTECTIVE EQUIPMENT

SKIN PROTECTION:

Skin contact with refrigerant may cause frostbite. General work clothing and gloves (leather) should provide adequate protection. If prolonged contact with liquid or gas is anticipated, insulated gloves constructed of PVA, neoprene or butyl rubber should be used. Any contaminated clothing should be promptly removed and washed before reuse.

EYE PROTECTION:

For normal conditions, wear safety glasses. Where there is reasonable probability of liquid contact, wear chemical safety goggles.

RESPIRATORY PROTECTION:

None generally required for adequately ventilated work situations. For accidental release or non-ventilated situations, or release into confined space, where the concentration may be above the PEL of 1,000 ppm, use a self-contained, NIOSH approved breathing apparatus or supplied air respirator. For escape: use the former or a NIOSH approved gas mask with organic vapor canister.

ADDITIONAL RECOMMENDATIONS:

Where contact with liquid is likely, such as in a spill or leak, impervious boots and clothing should be worn. High dose-level warning signs are recommended for areas of principle exposure. Provide eyewash stations and quick-drench shower facilities at convenient locations. For tank cleaning operations, see OSHA regulations, 29 CFR 1910.132 and 29 CFR 1910.133.

EXPOSURE GUIDELINES**INGREDIENT NAME**

1,1,1,2-Tetrafluoroethane

ACGIH TLV

None

OSHA PEL

None

OTHER LIMIT

*1000 ppm TWA (8hr)

* = Workplace Environmental Exposure Level (AIHA)

OTHER EXPOSURE LIMITS FOR POTENTIAL DECOMPOSITION PRODUCTS:

Hydrogen Fluoride: ACGIH TLV: 2 ppm ceiling, 0.5 ppm TLV-TWA

9. PHYSICAL AND CHEMICAL PROPERTIES

APPEARANCE:	Clear, colorless liquid and vapor	
PHYSICAL STATE:	Gas at ambient temperatures	
MOLECULAR WEIGHT:	102	
CHEMICAL FORMULA:	F_3CCH_2F	
ODOR:	Faint ethereal odor	
SPECIFIC GRAVITY (water = 1.0):	<1.22	
SOLUBILITY IN WATER (weight %):	0.15 wt%	
pH:	Neutral	
BOILING POINT:	-26.2°C (-15.1°F)	
FREEZING POINT:	-92.5°C (-141.9°F)	
VAPOR PRESSURE:	85.8 psia @ 70°F 213.4 psia @ 130°F	
VAPOR DENSITY (air = 1.0):	3.5	
EVAPORATION RATE:	>1	COMPARED TO: $CCl_4 = 1$
% VOLATILES:	100	
FLASH POINT:	Not applicable	

(Flash point method and additional flammability data are found in Section 5.)

10. STABILITY AND REACTIVITY

NORMALLY STABLE? (CONDITIONS TO AVOID):

The product is stable.

Do not mix with oxygen or air above atmospheric pressure. Any source of high temperatures, such as lighted cigarettes, flames, hot spots or welding may yield toxic and/or corrosive decomposition products.

INCOMPATIBILITIES:

(Under specific conditions: e.g. very high temperatures and/or appropriate pressures) – Freshly abraded aluminum surfaces (may cause strong exothermic reaction). Chemically reactive metals: potassium, calcium, powdered aluminum, magnesium, and zinc.

HAZARDOUS DECOMPOSITION PRODUCTS:

Halogens, halogen acids and possibly carbonyl halides.

HAZARDOUS POLYMERIZATION:

Will not occur.

11. TOXICOLOGICAL INFORMATION

IMMEDIATE (ACUTE) EFFECTS:

LC₅₀ : 4 hr. (rat) - > 500,000 ppm / Cardiac Sensitization threshold (dog) 80,000 ppm. NOEL – 50,000 ppm

DELAYED (SUBCHRONIC AND CHRONIC) EFFECTS:

Not mutagenic in four tests

Teratogenic NOEL (rat and rabbit) – 40,000 ppm

Subchronic inhalation (rat) NOEL – 50,000 ppm

Chronic NOEL – 10,000 ppm

OTHER DATA:

Metabolism <0.5% as CO₂ in tests at 50,000 ppm, late developing benign tumors were found.

12. ECOLOGICAL INFORMATION

Degradability (BOD):

R-134A is a gas at room temperature; therefore, it is unlikely to remain in water.

Octanol Water Partition Coefficient: Log P_{ow} = 1.06

13. DISPOSAL CONSIDERATIONS

RCRA

Is the unused product a RCRA hazardous waste if discarded?
If yes, the RCRA ID number is:

Not a hazardous waste
Not applicable

OTHER DISPOSAL CONSIDERATIONS:

Disposal must comply with federal, state, and local disposal or discharge laws. R-134A is subject to U.S. Environmental Protection Agency Clean Air Act Regulations Section 608 in 40 CFR Part 82 regarding refrigerant recycling.

The information offered here is for the product as shipped. Use and/or alterations to the product such as mixing with other materials may significantly change the characteristics of the material and alter the RCRA classification and the proper disposal method.

14. TRANSPORT INFORMATION

US DOT PROPER SHIPPING NAME: 1,1,1,2-Tetrafluoroethane
US DOT HAZARD CLASS: 2.2
US DOT PACKING GROUP: Not applicable
US DOT ID NUMBER: UN3159

For additional information on shipping regulations affecting this material, contact the information number found in Section 1.

15. REGULATORY INFORMATION**TOXIC SUBSTANCES CONTROL ACT (TSCA)**

TSCA INVENTORY STATUS: Listed on the TSCA inventory

OTHER TSCA ISSUES: None

SARA TITLE III / CERCLA

“Reportable Quantities” (RQs) and/or “Threshold Planning Quantities” (TPQs) exist for the following ingredients.

INGREDIENT NAME**SARA / CERCLA RQ (lb.)****SARA EHS TPQ (lb.)**

No ingredients listed in this section

Spills or releases resulting in the loss of any ingredient at or above its RQ requires immediate notification to the National Response Center [(800) 424-8802] and to your Local Emergency Planning Committee.

SECTION 311 HAZARD CLASS: IMMEDIATE
PRESSURE

SARA 313 TOXIC CHEMICALS:

The following ingredients are SARA 313 “Toxic Chemicals”. CAS numbers and weight percents are found in Section 2.

INGREDIENT NAME**COMMENT**

No ingredients listed in this section

STATE RIGHT-TO-KNOW

In addition to the ingredients found in Section 2, the following are listed for state right-to-know purposes.

INGREDIENT NAME**WEIGHT %****COMMENT**

No ingredients listed in this section

ADDITIONAL REGULATORY INFORMATION:

R-134A is subject to U.S. Environmental Protection Agency Clean Air Act Regulations at 40 CFR Part 82.

WARNING: DO NOT vent to the atmosphere. To comply with provisions of the U.S. Clean Air Act, any residual must be recovered. **Contains 1,1,1,2-Tetrafluoroethane (HFC-134a)**, a greenhouse gas which may contribute to global warming.

WHMIS CLASSIFICATION (CANADA):

This product has been evaluated in accordance with the hazard criteria of the CPR and the MSDS contains all the information required by the CPR.

FOREIGN INVENTORY STATUS:

Canada – Listed on DSL
EU - EINECS # 223770

16. OTHER INFORMATION

CURRENT ISSUE DATE: December, 2008

PREVIOUS ISSUE DATE: August, 2007

OTHER INFORMATION: HMIS Classification: Health – 1, Flammability – 1, Reactivity – 0
NFPA Classification: Health – 2, Flammability – 1, Reactivity – 0
ANSI/ASHRAE 34 Safety Group – A1
UL Classified

Regulatory Standards:

1. OSHA regulations for compressed gases: 29 CFR 1910.101
2. DOT classification per 49 CFR 172.101

Toxicity information per PAFT Testing

17. DISCLAIMER

National Refrigerants, Inc. believes that the information and recommendations contained herein (including data and statements are accurate as of the date hereof. NO WARRANTY OF FITNESS FOR ANY PARTICULAR PURPOSE, WARRANTY OF MERCHANTABILITY, OR ANY OTHER WARRANTY, EXPRESSED OR IMPLIED, IS MADE CONCERNING THE INFORMATION PROVIDED HEREIN. The information provided herein relates only to the specific product designated and may not be valid where such product is used in combination with any other methods of use of the product and of the information referred to herein are beyond the control of National Refrigerants. National Refrigerants expressly disclaims any and all liability as to any results obtained or arising from any use of the product or reliance on such information.

5c. Safety data sheet for R-404A

Material Safety Data Sheet

R-404A

1. CHEMICAL PRODUCT AND COMPANY IDENTIFICATION

PRODUCT NAME: R-404A
DISTRIBUTOR: National Refrigerants, Inc.
661 Kenyon Avenue
Bridgeton, New Jersey 08302

FOR MORE INFORMATION CALL:
(Monday-Friday, 8:00am-5:00pm)
1-800-262-0012

IN CASE OF EMERGENCY CALL:
CHEMTREC: 1-800-424-9300

2. COMPOSITION / INFORMATION ON INGREDIENTS

<u>INGREDIENT NAME</u>	<u>CAS NUMBER</u>	<u>WEIGHT %</u>
Pentafluoroethane (HFC-125)	354-33-6	44
1,1,1-Trifluoroethane (HFC-143a)	420-46-2	52
1,1,1,2-Tetrafluoroethane (HFC-134a)	811-97-2	4

Trace impurities and additional material names not listed above may also appear in Section 15 toward the end of the MSDS. These materials may be listed for local "Right-To-Know" compliance and for other reasons.

3. HAZARDS IDENTIFICATION

EMERGENCY OVERVIEW: Colorless, volatile liquid with ethereal and faint sweetish odor. Non-flammable material. Overexposure may cause dizziness and loss of concentration. At higher levels, CNS depression and cardiac arrhythmia may result from exposure. Vapors displace air and can cause asphyxiation in confined spaces. At higher temperatures, (>250°C), decomposition products may include Hydrofluoric Acid (HF) and carbonyl halides.

POTENTIAL HEALTH HAZARDS

SKIN: Irritation would result from a defatting action on tissue. Liquid contact could cause frostbite.

EYES: Liquid contact can cause severe irritation and frostbite. Mist may irritate.

INHALATION: R-404A is low in acute toxicity in animals. When oxygen levels in air are reduced to 12-14% by displacement, symptoms of asphyxiation, loss of coordination, increased pulse rate and deeper respiration will occur. At high levels, cardiac arrhythmia may occur.

INGESTION: Ingestion is unlikely because of the low boiling point of the material. Should it occur, discomfort in the gastrointestinal tract from rapid evaporation of the material and consequent evolution of gas would result. Some effects of inhalation and skin exposure would be expected.

DELAYED EFFECTS: None known.

Ingredients found on one of the OSHA designated carcinogen lists are listed below.

<u>INGREDIENT NAME</u>	<u>NTP STATUS</u>	<u>IARC STATUS</u>	<u>OSHA LIST</u>
No ingredients listed in this section			

4. FIRST AID MEASURES

SKIN: Promptly flush skin with water until all chemical is removed. If there is evidence of frostbite, bathe (do not rub) with lukewarm (not hot) water. If water is not available, cover with a clean, soft cloth or similar covering. Get medical attention if symptoms persist.

EYES: Immediately flush eyes with large amounts of water for at least 15 minutes (in case of frostbite water should be lukewarm, not hot) lifting eyelids occasionally to facilitate irrigation. Get medical attention if symptoms persist.

INHALATION: Immediately remove to fresh air. If breathing has stopped, give artificial respiration. Use oxygen as required, provided a qualified operator is available. Get medical attention. Do not give epinephrine (adrenaline).

INGESTION: Ingestion is unlikely because of the physical properties and is not expected to be hazardous. Do not induce vomiting unless instructed to do so by a physician.

ADVICE TO PHYSICIAN: Because of the possible disturbances of cardiac rhythm, catecholamine drugs, such as epinephrine, should be used with special caution and only in situations of emergency life support. Treatment of overexposure should be directed at the control of symptoms and the clinical conditions.

5. FIRE FIGHTING MEASURES

FLAMMABLE PROPERTIES

FLASH POINT:	Gas, not applicable per DOT regulations
FLASH POINT METHOD:	Not applicable
AUTOIGNITION TEMPERATURE:	<750°C
UPPER FLAME LIMIT (volume % in air):	None*
LOWER FLAME LIMIT (volume % in air):	None*
	*Based on ASHRAE Standard 34 with match ignition
FLAME PROPAGATION RATE (solids):	Not applicable
OSHA FLAMMABILITY CLASS:	Not applicable

EXTINGUISHING MEDIA:

Use any standard agent – choose the one most appropriate for type of surrounding fire (material itself is not flammable)

UNUSUAL FIRE AND EXPLOSION HAZARDS:

R-404A is not flammable at ambient temperatures and atmospheric pressure. However, this material will become combustible when mixed with air under pressure and exposed to strong ignition sources. Contact with certain reactive metals may result in formation of explosive or exothermic reactions under specific conditions (e.g. very high temperatures and/or appropriate pressures).

SPECIAL FIRE FIGHTING PRECAUTIONS/INSTRUCTIONS:

Firefighters should wear self-contained, NIOSH-approved breathing apparatus for protection against possible toxic decomposition products. Proper eye and skin protection should be provided. Use water spray to keep fire-exposed containers cool.

6. ACCIDENTAL RELEASE MEASURES

IN CASE OF SPILL OR OTHER RELEASE:

(Always wear recommended personal protective equipment.)

Evacuate unprotected personnel. Protected personnel should remove ignition sources and shut off leak, if without risk, and provide ventilation. Unprotected personnel should not return until air has been tested and determined safe, including low-lying areas.

Spills and releases may have to be reported to Federal and/or local authorities. See Section 15 regarding reporting requirements.

7. HANDLING AND STORAGE

NORMAL HANDLING:

(Always wear recommended personal protective equipment.)

Avoid breathing vapors and liquid contact with eyes, skin or clothing. Do not puncture or drop cylinders, expose them to open flame or excessive heat. Use authorized cylinders only. Follow standard safety precautions for handling and use of compressed gas cylinders.

R-404A should not be mixed with air above atmospheric pressure for leak testing or any other purpose.

STORAGE RECOMMENDATIONS:

Store in a cool, well-ventilated area of low fire risk and out of direct sunlight. Protect cylinder and its fittings from physical damage. Storage in subsurface locations should be avoided. Close valve tightly after use and when empty.

8. EXPOSURE CONTROLS / PERSONAL PROTECTION

ENGINEERING CONTROLS:

Provide local ventilation at filling zones and areas where leakage is probable. Mechanical (general) ventilation may be adequate for other operating and storage areas.

PERSONAL PROTECTIVE EQUIPMENT

SKIN PROTECTION:

Skin contact with refrigerant may cause frostbite. General work clothing and gloves (leather) should provide adequate protection. If prolonged contact with the liquid or gas is anticipated, insulated gloves constructed of PVA, neoprene or butyl rubber should be used. Any contaminated clothing should be promptly removed and washed before reuse.

EYE PROTECTION:

For normal conditions, wear safety glasses. Where there is reasonable probability of liquid contact, wear chemical safety goggles.

RESPIRATORY PROTECTION:

None generally required for adequately ventilated work situations. For accidental release or non-ventilated situations, or release into confined space, where the concentration may be above the PEL of 1,000 ppm, use a self-contained, NIOSH-approved breathing apparatus or supplied air respirator. For escape: use the former or a NIOSH-approved gas mask with organic vapor canister.

ADDITIONAL RECOMMENDATIONS:

Where contact with liquid is likely, such as in a spill or leak, impervious boots and clothing should be worn. High dose-level warning signs are recommended for areas of principle exposure. Provide eyewash stations and quick-drench shower facilities at convenient locations. For tank cleaning operations, see OSHA regulations, 29 CFR 1910.132 and 29 CFR 1910.133.

EXPOSURE GUIDELINES

<u>INGREDIENT NAME</u>	<u>ACGIH TLV</u>	<u>OSHA PEL</u>	<u>OTHER LIMIT</u>
Pentafluoroethane	None	None	*1000 ppm TWA (8hr)
1,1,1-Trifluoroethane	None	None	*1000 ppm TWA (8hr)
1,1,1,2-Tetrafluoroethane	None	None	*1000 ppm TWA (8hr)

* = Workplace Environmental Exposure Level (AIHA)

OTHER EXPOSURE LIMITS FOR POTENTIAL DECOMPOSITION PRODUCTS:

Hydrogen Fluoride: ACGIH TLV: 3 ppm ceiling

9. PHYSICAL AND CHEMICAL PROPERTIES

APPEARANCE:	Clear, colorless liquid and vapor	
PHYSICAL STATE:	Gas at ambient temperatures	
MOLECULAR WEIGHT:	120	
CHEMICAL FORMULA:	CHF ₂ CF ₃ , CH ₃ CF ₃ , CH ₂ FCF ₃	
ODOR:	Faint ethereal odor	
SPECIFIC GRAVITY (water = 1.0):	1.08 @ 21.1°C (70°F)	
SOLUBILITY IN WATER (weight %):	Unknown	
pH:	Neutral	
BOILING POINT:	-47.8°C (-54.0°F)	
FREEZING POINT:	Not Determined	
VAPOR PRESSURE:	182.9 psia @ 70°F 370.9 psia @ 130°F	
VAPOR DENSITY (air = 1.0):	3.43	
EVAPORATION RATE:	>1	COMPARED TO: CC1 ₄ = 1
% VOLATILES:	100	
FLASH POINT:	Not applicable	
(Flash point method and additional flammability data are found in Section 5.)		

10. STABILITY AND REACTIVITY

NORMALLY STABLE? (CONDITIONS TO AVOID):

The product is stable.

Do not mix with oxygen or air above atmospheric pressure. Any source of high temperature, such as lighted cigarettes, flames, hot spots or welding may yield toxic and/or corrosive decomposition products.

INCOMPATIBILITIES:

(Under specific conditions: e.g. very high temperatures and/or appropriate pressures) – Freshly abraded aluminum surfaces (may cause strong exothermic reaction). Chemically active metals: potassium, calcium, powdered aluminum, magnesium and zinc.

HAZARDOUS DECOMPOSITION PRODUCTS:

Halogens, halogen acids and possibly carbonyl halides.

HAZARDOUS POLYMERIZATION:

Will not occur.

11. TOXICOLOGICAL INFORMATION

IMMEDIATE (ACUTE) EFFECTS:

HFC-125: LC₅₀ : 4 hr. (rat) - > 800,000 ppm / Cardiac Sensitization threshold (dog) 75,000 ppm

HFC-143a: LC₅₀ : 4hr. (rat) - > 540,000 ppm / Cardiac Sensitization threshold (dog) > 250,000 ppm

HFC-134a: LC₅₀ : 4hr. (rat) - > 500,000 ppm / Cardiac Sensitization threshold (dog) > 80,000 ppm

DELAYED (SUBCHRONIC AND CHRONIC) EFFECTS:

HFC-125: Teratogenic NOEL (rat and rabbit) – 50,000 ppm

Subchronic inhalation (rat) NOEL - $\geq$ 50,000 ppm / Chronic NOEL – 10,000 ppm

HFC-143a: Teratogenic NOEL (rat and rabbit) – 50,000 ppm

Subchronic inhalation (rat) NOEL - $\geq$ 50,000 ppm

HFC-134a: Teratogenic NOEL (rat and rabbit) – 40,000 ppm

Subchronic inhalation (rat) NOEL – 50,000 ppm / Chronic NOEL – 10,000 ppm

OTHER DATA:

HFC-125, HFC-134a: Not active in four genetic studies

HFC-143a: Not active in two genetic studies

12. ECOLOGICAL INFORMATION

Degradability (BOD):

R-404A is a gas at room temperature; therefore, it is unlikely to remain in water.

Octanol Water Partition Coefficient: Unknown for mixture

13. DISPOSAL CONSIDERATIONS

RCRA

Is the unused product a RCRA hazardous waste if discarded?
If yes, the RCRA ID number is:

Not a hazardous waste.
Not applicable.

OTHER DISPOSAL CONSIDERATIONS:

Disposal must comply with federal, state, and local disposal or discharge laws. R-404A is subject to U.S. Environmental Protection Agency Clean Air Act Regulations Section 608 in 40 CFR Part 82 regarding refrigerant recycling.

The information offered here is for the product as shipped. Use and/or alterations to the product such as mixing with other materials may significantly change the characteristics of the material and alter the RCRA classification and the proper disposal method.

14. TRANSPORT INFORMATION

US DOT PROPER SHIPPING NAME: Refrigerant gas R 404A
US DOT HAZARD CLASS: 2.2
US DOT PACKING GROUP: Not applicable
US DOT ID NUMBER: UN3337

For additional information on shipping regulations affecting this material, contact the information number found in Section 1.

15. REGULATORY INFORMATION

TOXIC SUBSTANCES CONTROL ACT (TSCA)

TSCA INVENTORY STATUS: Components listed on the TSCA inventory
OTHER TSCA ISSUES: None

SARA TITLE III / CERCLA

“Reportable Quantities” (RQs) and/or “Threshold Planning Quantities” (TPQs) exist for the following ingredients.

INGREDIENT NAME**SARA / CERCLA RQ (lb.)****SARA EHS TPO (lb.)**

No ingredients listed in this section

Spills or releases resulting in the loss of any ingredient at or above its RQ requires immediate notification to the National Response Center [(800) 424-8802] and to your Local Emergency Planning Committee.

SECTION 311 HAZARD CLASS: IMMEDIATE PRESSURE

SARA 313 TOXIC CHEMICALS:

The following ingredients are SARA 313 “Toxic Chemicals”. CAS numbers and weight percents are found in Section 2.

INGREDIENT NAME**COMMENT**

No ingredients listed in this section

STATE RIGHT-TO-KNOW

In addition to the ingredients found in Section 2, the following are listed for state right-to-know purposes.

INGREDIENT NAME**WEIGHT %****COMMENT**

No ingredients listed in this section

ADDITIONAL REGULATORY INFORMATION:

R-404A is subject to U.S. Environmental Protection Agency Clean Air Act Regulations at 40 CFR Part 82.

WARNING: Contains Pentafluoroethane (HFC-125), 1,1,1-trifluoroethane, tetrafluoroethane, greenhouse gases which may contribute to global warming. **Do not vent** to the atmosphere. To comply with provisions of the U.S. Clean Air Act, any residual must be recovered.

WHMIS CLASSIFICATION (CANADA):

This product has been evaluated in accordance with the hazard criteria of the CPR and the MSDS contains all the information required by the CPR.

FOREIGN INVENTORY STATUS:

EU – EINECS # 2065578 – HFC-125
2069965 – HFC-143a
223770 – HFC-134a

16. OTHER INFORMATION

CURRENT ISSUE DATE: December, 2008

PREVIOUS ISSUE DATE: August, 2007

OTHER INFORMATION: HMIS Classification: Health – 1, Flammability – 1, Reactivity – 0
NFPA Classification: Health – 2, Flammability – 1, Reactivity – 0
ANSI / ASHRAE 34 Safety Group – A1

Regulatory Standards:

1. OSHA regulations for compressed gases: 29 CFR 1910.101
2. DOT classification per 49 CFR 172.101

Toxicity information per PAFT Testing



17. DISCLAIMER

National Refrigerants, Inc. believes that the information and recommendations contained herein (including data and statements) are accurate as of the date hereof. NO WARRANTY OF FITNESS FOR ANY PARTICULAR PURPOSE, WARRANTY OF MERCHANTABILITY, OR ANY OTHER WARRANTY, EXPRESSED OR IMPLIED, IS MADE CONCERNING THE INFORMATION PROVIDED HEREIN. The information provided herein relates only to the specific product designated and may not be valid where such product is used in combination with any other methods or use of the product and of the information referred to herein are beyond the control of National Refrigerants. National Refrigerants expressly disclaims any and all liability as to any results obtained or arising from any use of the product or reliance on such information.

MSDS: R-404A
Current Issue Date: December, 2008

Page 7 of 7

5e. Controller G-214

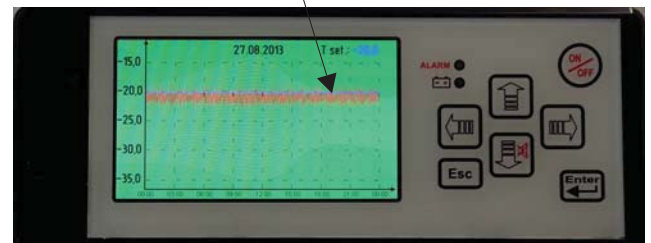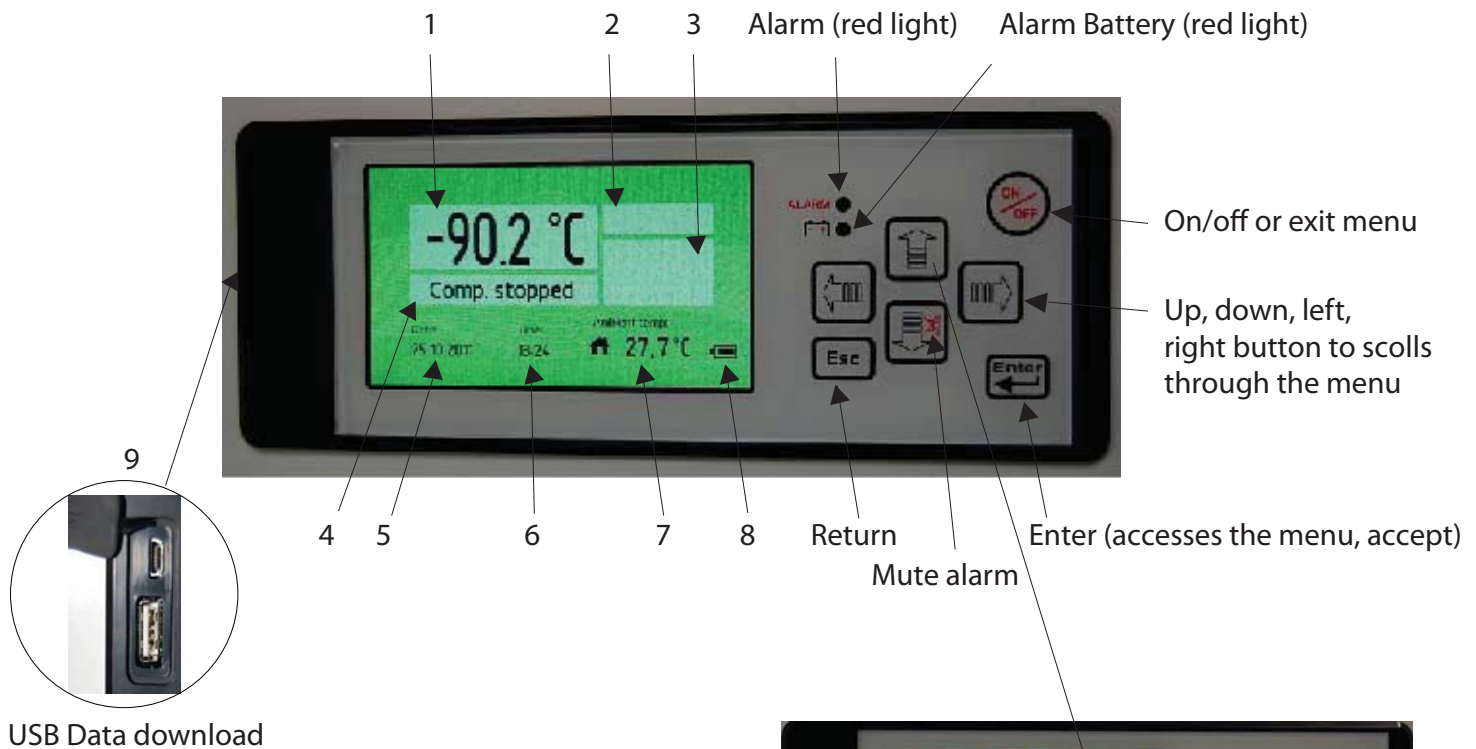

Control Panel Description


1. Temperature indicator
2. Alarm indicator
3. Alarm-icons
4. Compressor indication
5. Date indication
6. Hour indication
7. Ambient temperature indicator
8. Battery level
9. USB-Connection

Push the **UP** button in order to generate a graph.


Push the **LEFT** button in order to return up to 10 days back.

Push the **Esc** button in order to exit the menu.


Setting icons

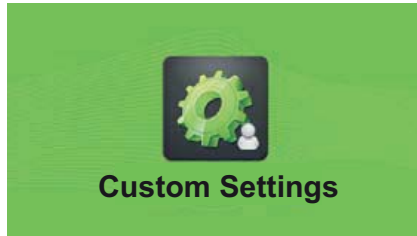


Custom Settings

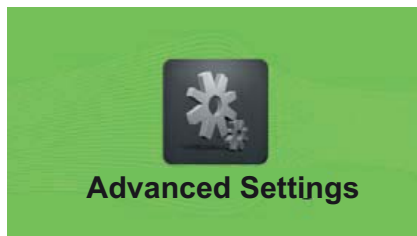


Advanced Settings

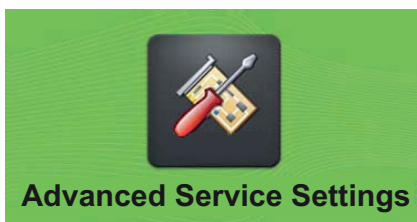


Advanced Service Settings

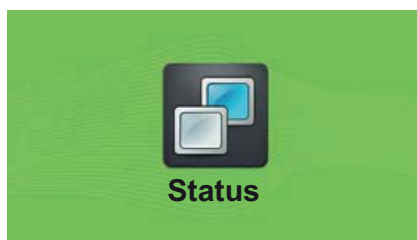

Status

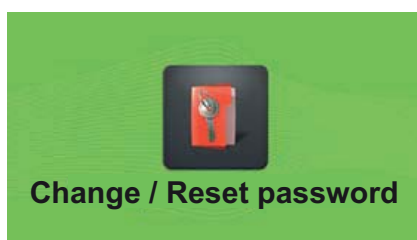

Change / Reset password

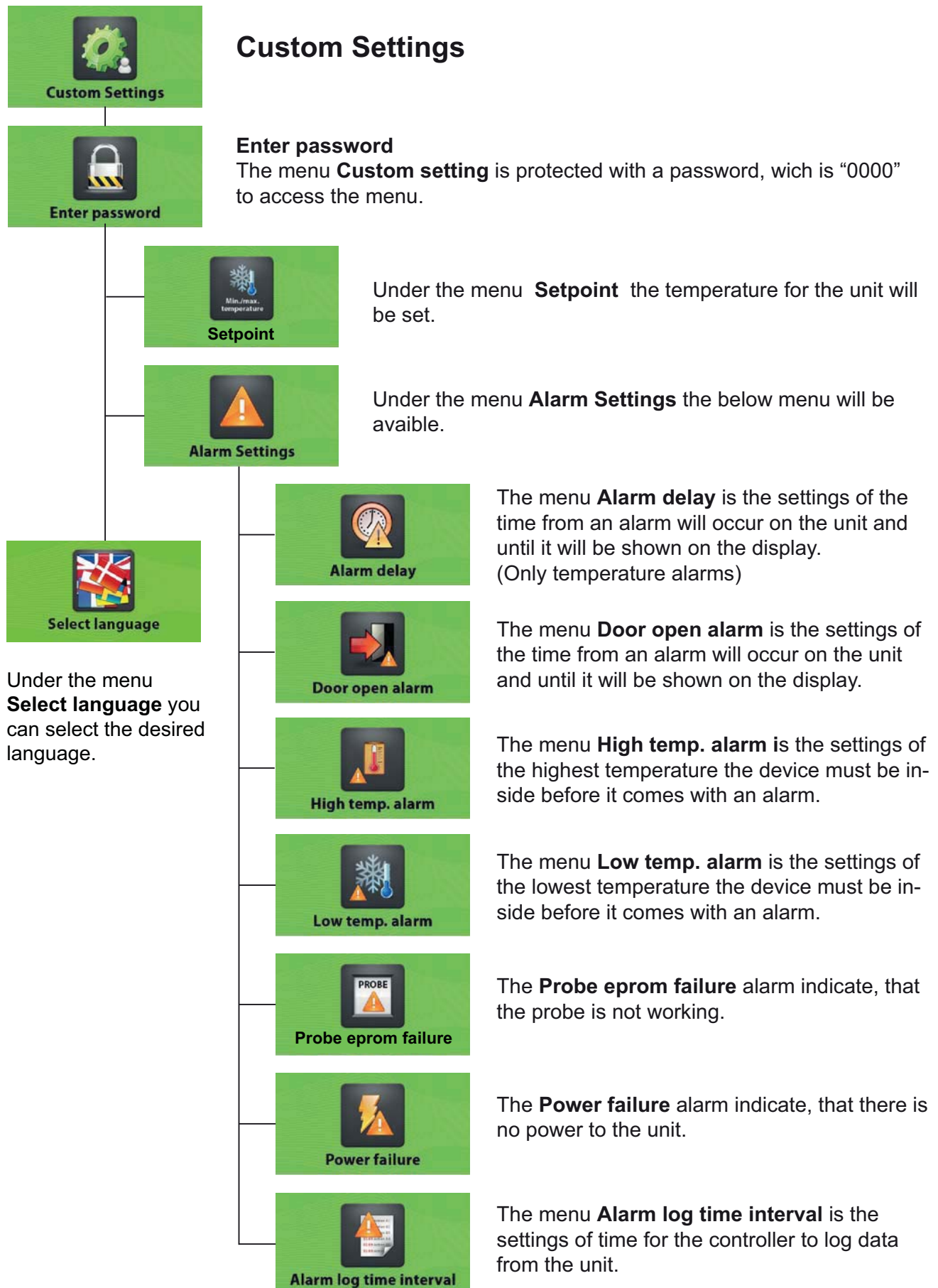




Advanced Settings



Enter password

The menu **Advanced Settings** is protected with a password which is "0000" to access the menu.

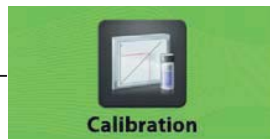


The menu **Calibration** is the settings for offset of the temperature in the display.



The menu **Automatic defrost** is the setting of the time where a defrost automatic will be performed on the unit.




Under the menu **Manual defrost** the defrost can be started immediately.

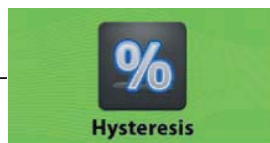


The menu **Hysteresis** is the setting of the difference from the setpoint where the compressor will start and stop.



The menu **Temp. range limits** is the settings of the maximum and minimum temperature the unit can be used for.

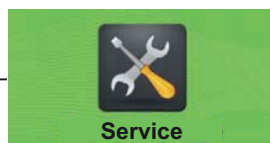

→ **Service**

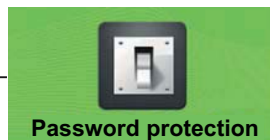

The menu **Password protection** is the setting of the protection for the unit not to be turned off and is protected with a password.



The menu **Set date / time** is the setting of date and time.



The menu **Ambient temp. settings** is the adjustment of the ambient temperature.



Service

Service Settings

Under the menu **Service**, there is information about the unit, which is important for service on the unit.



Compressor hours

The menu **Compressor hours** shows, how many hours the compressor has run.



Fan hours

The menu **Fan hours** shows, how many hours the fan has run.



Probe

Under the menu **Probe** the temperature for the sensor mounted in the unit is displayed.



Probe inside unit

Under the menu **Probe inside unit** the temperature of the sensor mounted inside the unit is displayed. There is mounted a sensor as standard in the unit, but it is possible to mount 3 sensors inside the device. (Controller version)



Probe Compressor

Under the item **Probe Compressor** the temperature of the sensor mounted on the compressor is displayed. (Not all models)



Probe evaporator

Under the point **Probe evaporator** the temperature of the sensor, that is located on the condenser, is displayed.



Advanced Service Settings

Advanced Service Settings

This menu is only used by the service department at the manufacturer.



Status

Status



Current alarms

The **Current alarms** shows the alarm, which have been on the unit with data and time.



Temp. last 24 hours

The **Temp. last 24 hours** shows the temperature of the unit the last 24 hours.



Change/Reset password

Change / Reset password



Change password

Under **Change password** it is possible to change the password for **Customer Settings**, **Advanced Settings** and **Advanced Service Settings**.



Reset password

Under **Reset password** it is possible to reset the password for the **Customer Settings**, **Advanced Settings** and **Advanced Service Settings**.

Contact the manufacturer to get the password for reset password.

Download Data



Place the USB pen drive in the connection on the display and press the button right for **read data** from the unit to the USB pen.



The **Read data** will be shown in the display when the data is loading to the USB pen drive.



The **Transfer complete** will be shown on the display when the loading is finished and the USB pen can be removed from the display.

Open the 2 files:
data00 and **param00**
on a computer in Excel or similar.


Upload Data (Only for new settings from the manufactory)



Turn-off the display before start.
Place the USB pen drive in the connection on the display and press the button left for upload parameters from the USB pen to the unit.



The **Data upload** will be shown in the display when the data is loading to the unit.

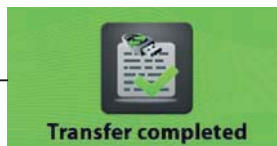

The **Transfer complete** will be shown on the display when the loading is finished and the USB pen can be removed from the display.

Alarms on display




Alarm icons




The **Door open alarm** indicate, that the door is not correct closed



The **High temperature alarm** indicate, that the temperature inside the unit is higher than permitted from setting of the unit.

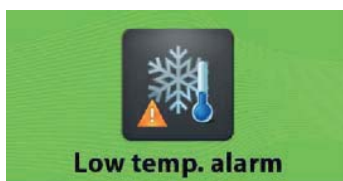

The **Low temperature alarm** indicate, that the temperature inside the unit is lower than permitted from setting of the unit.

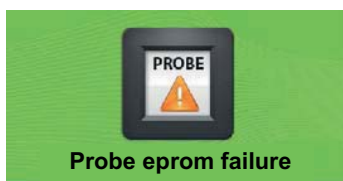

The **Probe eprom failure** alarm indicate, that the probe is not working.



The **Power failure** alarm indicate, that there is no power to the unit.



The **Low battery** alarm indicate, that the power on the battery is to low.

5e. Default settings

Description	Arctiko default settings	Customer settings
Custom Settings (Password)	0000	
Set point		
Refrigerator	+4,0	
Alarm Settings:		
Alarm delay		
Refrigerator	15 min	
Door open alarm		
Refrigerator	Enable	
High temp. alarm		
Refrigerator	+6	
Low temp. alarm		
Refrigerator	+2	
Probe/eprom failure		
Refrigerator	Enable	
Power failure		
Refrigerator	Enable	
Alarm log time interval		
Refrigerator	1	

Description	Arctiko default settings	Customer settings
Advanced Settings (Password)	0000	
Calibration		
Refrigerator	0,0	
Automatic defrost		
Refrigerator	06	
Manual defrost		
Refrigerator	OK	
Hysteresis		
Refrigerator	1,0	
Temp. range limits		
Refrigerator Max.	+10,0	
Refrigerator Min.	+1,0	

ARCTIKO

The art of simplicity®



www.arctiko.com

Lammefjordsvej 5
DK-6715 Esbjerg N
Denmark
Tel. +45 7020 0328
Fax. +45 7020 0329
info@arctiko.com
www.arctiko.com